US_EN

motorola


MODEL:
MBP69SN

Quick Start Guide

For a full explanation of all features and instructions, please refer to the User's Guide (available for download from www.motorolastore.com).


1. Setting Up the Smart Ear Thermometer



- The key to obtaining accurate readings is a clean sensor.
To clean the sensor, use a cotton swab slightly moistened with alcohol to gently wipe its surface, and then immediately wipe dry using a clean cotton swab.
Wait at least 20 minutes to make sure that the sensor is completely dry before taking temperature.
- Make sure the ear canal is clean. If you have cleaned the ear, please wait for 5 to 10 minutes before measuring temperature.
- If the temperature of the storage area differs greatly from that of the measuring area, please wait until the thermometer temperature has equalized to the room temperature, about 30 minutes, prior to use.
- Please open the battery cover and install 2 AAA batteries, which can be found inside the packaging. Once batteries are installed, the thermometer will then turn on automatically and the LCD display will appear as in the figure to the right. Press any button to turn the device back on when it turns itself off.


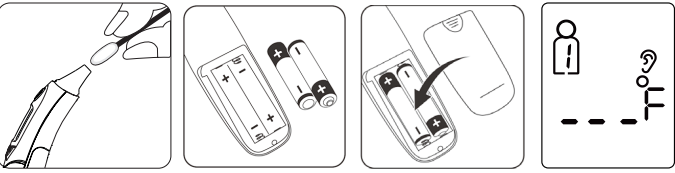
HOW TO MEASURE EAR TEMPERATURE

- Make sure the probe is clean before measuring. Remove the protective cap prior to use. Place the thermometer in your hand with your thumb or forefinger on the [] button.
 - Do not press the button yet.
 - Make sure the probe is clean before measuring.

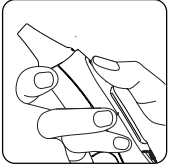- Insert the probe very gently and slowly into the ear canal.
 - Hold the thermometer such that the probe is oriented straight in the direction of the eardrum.
 - It is important to point the thermometer's probe window toward the eardrum properly inside the ear.
 - Hold the outer ear and gently pull it toward the rear of the head to straighten the ear canal.
 - Always insert the thermometer into the same ear in the same direction and at the same depth.

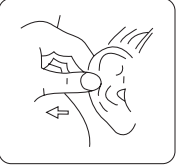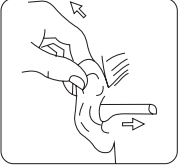- Press the [] button to start measurement until you hear two "DIDI" beep. (approx. 1 second). The measurement is complete.
- Remove the thermometer from the ear canal. The LCD displays the measured temperature. You must wait at least 6 seconds or until the green or red LED turns off before another measurement can be taken. When the cleaning symbol flashes on the LED screen, please clean the infra-red sensor. Press any button to stop flashing.
- Press the [] button for 2 seconds to mute the device.
- If no button is pressed the unit will switch off after 1 minute.

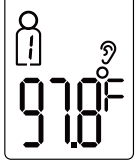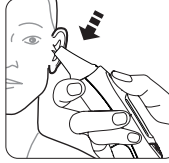Please take note of the following minimum system requirements:
Smartphones/ Tablets: Android™ and iOS® only.
Minimum requirement: iOS® 8.0, Android™ 4.4.2 or above.
For more information, please visit: <http://hubblehome.com/requirements>

Install Hubble Connect for Smart Nursery App on your Android™ Devices/ iPhone®/ iPad®.

Download **Hubble Connect for Smart Nursery** App from the App Store for iOS devices or from the Google Play™ Store for Android™ devices.

Pair the Smart Ear Thermometer with your Android™ Device/ iPhone®/ iPad®.

- Run the Hubble Connect for Smart Nursery app on your Android™ Device/iPhone®/iPad®.
- Turn on Bluetooth on your Android™ Device iPhone®/ iPad®.
- Follow the in-app instructions to create your Smart Nursery account and connect to the Thermometer via Bluetooth.
- Log in to your Hubble Connect for Smart Nursery account, and tap on "Smart Ear Thermometer" to enter the set-up mode.
- Tap on "Skip to Setup", then tap "Next". Make sure your thermometer is switched on, then tap next to search for the thermometer
- Search for the device, and when the device is found, tap "MBP69-XXX".
- It will take a few seconds for the thermometer to connect to your smart device. Tap "Continue".
- Tap "Get Started" when the thermometer is successfully connected to the Hubble Connect for Smart Nursery app.

Note
Should you fail to complete the set-up process, repeat the above steps (5).

2. General Information

If your product is not working properly, read this Quick Start Guide or the User's Guide.
Contact Customer Service:
US and Canada, 1-888-331-3383; website: www.motorolastore.com

Consumer Products and Accessories Limited Warranty ("Warranty")
Thank you for purchasing this Motorola branded product manufactured under license by Binatone Electronics International Ltd. ("BINATONE").

What Does This Warranty Cover?
Subject to the exclusions contained below, BINATONE warrants that this Motorola branded product ("Product") or certified accessory ("Accessory") sold for use with this product is manufactured to be free from defects in materials and workmanship under normal consumer usage for the period outlined below. This Warranty is your exclusive warranty and is not transferable.

Who is covered?
This Warranty extends only to the first consumer purchaser, and is not transferable.

What Will BINATONE Do?
BINATONE or its authorized distributor, at its option and within a commercially reasonable time, will at no charge repair or replace any Products or Accessories that does not conform to this Warranty. We may use functionally equivalent reconditioned/refurbished/pre-owned or new Products, Accessories or parts.

What Other Limitations Are There?

ANY IMPLIED WARRANTIES, INCLUDING WITHOUT LIMITATION THE IMPLIED WARRANTIES OF MERCHANTABILITY AND FITNESS FOR A PARTICULAR PURPOSE, SHALL BE LIMITED TO THE DURATION OF THIS LIMITED WARRANTY, OTHERWISE THE REPAIR OR REPLACEMENT PROVIDED UNDER THIS EXPRESS LIMITED WARRANTY IS THE EXCLUSIVE REMEDY OF THE CONSUMER, AND IS PROVIDED IN LIEU OF ALL OTHER WARRANTIES, EXPRESS OR IMPLIED.

IN NO EVENT SHALL MOTOROLA OR BINATONE BE LIABLE, WHETHER IN CONTRACT OR TORT (INCLUDING NEGLIGENCE) FOR DAMAGES IN EXCESS OF THE PURCHASE PRICE OF THE PRODUCT OR ACCESSORY, OR FOR ANY INDIRECT, INCIDENTAL, SPECIAL OR CONSEQUENTIAL DAMAGES OF ANY KIND, OR LOSS OF REVENUE OR PROFITS, LOSS OF BUSINESS, LOSS OF INFORMATION OR OTHER FINANCIAL LOSS ARISING OUT OF OR IN CONNECTION WITH THE ABILITY OR INABILITY TO USE THE PRODUCTS OR ACCESSORIES TO THE FULL EXTENT THESE DAMAGES MAY BE DISCLAIMED BY LAW.

Some jurisdictions do not allow the limitation or exclusion of incidental or consequential damages, or limitation on the length of an implied warranty, so the above limitations or exclusions may not apply to you. This Warranty gives you specific legal rights, and you may also have other rights that vary from one jurisdiction to another.

Products Covered	Length of Coverage
Consumer Products	One (1) year from the date of the products, original purchase by the first consumer purchaser of the product.
Consumer Accessories	Ninety (90) days from the date of the Accessories, original purchase by the first consumer purchaser of the product.
Consumer Products and Accessories that are Repaired or Replaced	The balance of the original warranty or for Ninety (90) days from the date returned to the consumer, whichever is longer.

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Exclusions

Normal Wear and Tear. Periodic maintenance, repair and replacement of parts due to normal wear and tear are excluded from coverage.

Batteries. Only batteries whose fully charged capacity falls below 80% of their rated capacity and batteries that leak are covered by this Warranty.

Abuse & Misuse. Defects or damage that result from: (a) improper operation, storage, misuse or abuse, accident or neglect, such as physical damage (cracks, scratches, etc.) to the surface of the product resulting from misuse; (b) contact with liquid, water, rain, extreme humidity or heavy perspiration, sand, dirt or the like, extreme heat, or food; (c) use of the Products or Accessories for commercial purposes or subjecting the Product or Accessory to abnormal usage or conditions; or (d) other acts which are not the fault of MOTOROLA or BINATONE, are excluded from coverage.

Use of Non-Motorola branded Products and Accessories. Defects or damage that result from the use of Non-Motorola branded or certified Products or Accessories or other peripheral equipment are excluded from coverage.

Unauthorized Service or Modification. Defects or damages resulting from service, testing, adjustment, installation, maintenance, alteration, or modification in any way by someone other than MOTOROLA, BINATONE or its authorized service centers, are excluded from coverage.

Altered Products. Products or Accessories with (a) serial numbers or date tags that have been removed, altered or obliterated; (b) broken seals or that show evidence of tampering; (c) mismatched board serial numbers; or (d) non-conforming or non-Motorola branded housings, or parts, are excluded from coverage.

Communication Services. Defects, damages, or the failure of Products or Accessories due to any communication service or signal you may subscribe to or use with the Products or Accessories is excluded from coverage.

How to Obtain Warranty Service or Other Information ?

To obtain service or information, please call:

US and Canada, 1-888-331-3383; website: www.motorolastore.com

You will receive instructions on how to ship the Products or Accessories at your expense and risk, to a BINATONE Authorised Repair Center. To obtain service, you must include: (a) the Product or Accessory; (b) the original proof of purchase (receipt) which includes the date, place and seller of the Product; (c) if a warranty card was included in your box, a completed warranty card showing the serial number of the Product; (d) a written description of the problem; and, most importantly; (e) your address and telephone number.

These terms and conditions constitute the complete warranty agreement between you and BINATONE regarding the Products or Accessories purchased by you, and supersede any prior agreement or representations, including representations made in any literature publications or promotional materials issued by BINATONE or representations made by any agent, employee or staff of BINATONE, that may have been made in connection with the said purchase.

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Disposal of the Device (Environment)

At the end of the product life cycle, you should not dispose of this product with normal household waste. Take this product to a collection point for the recycling of electrical and electronic equipment. The symbol on the product, user's guide and/or box indicates this obligation.

By taking end-of-life products to a recycling point, you are ensuring that some of the product materials can be re-used. By reusing some parts or raw materials from used products, you make an important contribution to the protection of the environment.

Please contact your local authorities should you require additional information on the collection points in your area.

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IMPORTANT INFORMATION REGARDING ELECTROMAGNETIC COMPATIBILITY(EMC)

With the increased number of electronic devices such as PC's and mobile (cellular) telephones, medical devices in use may be susceptible to electromagnetic interference from other devices. Electromagnetic interference may result in incorrect operation of the medical device and create a potentially unsafe situation.

Medical devices should also not interfere with other devices.

In order to regulate the requirements for EMC (Electro Magnetic Compatibility) with the aim to prevent unsafe product situations, the IEC 60601-1-2:2014 standard has been implemented. This standard defines the levels of immunity to electromagnetic interferences as well as maximum levels of electromagnetic emissions for medical devices. Medical devices manufactured by Shenzhen Dongdixin Technology Co., Ltd. conform to this IEC 60601-1-2:2014 standard for both immunity and emissions. Nevertheless, special precautions need to be observed:

The use of accessories and cables other than those specified by manufacturer, with the exception of cables sold by manufacturer as replacement parts for internal components, may result in increased emission or decreased immunity of the device.

- The medical devices should not be used adjacent to or stacked with other equipment. In case adjacent or stacked use is necessary, the medical device should be observed to verify normal operation in the configuration in which it will be used.
- Refer to further guidance below regarding the EMC environment in which the device should be used.

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FCC and IC regulations (FCC Part 15)

Warning: Changes or modifications to this equipment not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

NOTE: This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to Part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications.

However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

Reorient or relocate the receiving antenna.
Increase the separation between the equipment and receiver.
Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
Consult the dealer or an experienced radio/TV technician for help.

This device complies with part 15 of the FCC Rules. Operation is subject to the following two conditions:
(1) This device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

Industry Canada. This device complies with the Industry Canada license-exempt RSS standard(s):
Operation is subject to the following Two conditions:(1) this device may not cause interference, and (2)

this device must accept any interference, including interference that may cause undesired operation of the device.

Download User Guide from website: www.motorolastore.com

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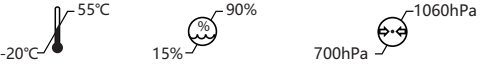
Manufactured, distributed or sold by Binatone Electronics International LTD., official licensee for this product. MOTOROLA and the Stylized M Logo are trademarks or registered trademarks of Motorola Trademark Holdings, LLC. and are used under license. Apple logo is a trademark of Apple Inc., registered in the U.S. and other countries. App Store is a service mark of Apple Inc. Google Play and Android are trademarks of Google Inc. Wi-Fi is a trademark of the Wi-Fi Alliance. All other trademarks are the property of their respective owners. © 2017 Motorola Mobility LLC. All rights reserved.

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