



breo Temple and Eye Massager
EM005

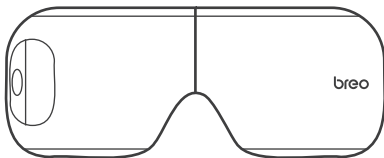

User Manual

breo[®] Temple and Eye Massager

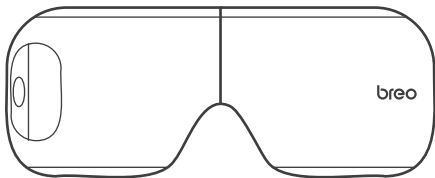

MENU

Attention	3
Indications for use	5
Product Structure	5
What's Included	6
How to Use	7
How to Connect Breo App	8
How to Wear	21
Charging	22
Caution	23
Symbol description	24
Cleaning and Maintenance	25
Troubleshooting	26
Guarantee	27
Product parameters	27
FCC Regulatory Compliance	28
Disposal	29
Customer Service	30



Thank you for your purchasing Breo Temple and Eye Massager!
Hope it helps you relax your eye muscle.
Please read this user's manual carefully prior to use and,
keep it properly for future reference.

Attention

This appliance can be used by children aged from 8 years and above and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved.

Children shall not play with the appliance.

Cleaning and user maintenance shall not be made by children without supervision.

Please read and comprehend all the contents of this manual and operate the product according to the instructions.

1. DO NOT use if you have had an eye operation, or if you have any eye diseases such as glaucoma, cataract, and amotio retinae.
2. DO NOT clean the product with water; clean it with cotton soaked in alcohol.
3. The appliance has a heated surface. Persons insensitive to heat must be careful when using the appliance.
4. Children, Juveniles, disabled people, and the elderly should only use this product when properly supervised.
5. The product is intended for personal use; it is not intended to replace medical care and medical devices.
6. DO NOT use if you have any swelling in the eye area.
7. DO NOT disassemble the product, DO NOT impact or squeeze its interior lithium battery or expose it to fire.
8. DO NOT use if the battery has been immersed in water.
9. DO NOT use in high temperatures or humid environments. Avoid direct exposure to the sun.
10. Consult your doctor if you have sensitive skin, allergic dermatitis, or orthopedic biomaterials.
11. Remove contact lenses, hairpins, and hair accessories before using.
12. Start the eye massager after putting it on in order to avoid damage to the air bag due to over pressure.
13. It must only be supplied at safety extra low voltage corresponding to the marking on the appliance.
14. You must have your eyes closed while massaging.
15. DO NOT use sharp objects to scratch this product
16. This appliance contains batteries that are non-replaceable.


Indications for use:

Breo Temple and Eye Massager is intended to relax eye muscle and relieve minor eyes muscle aches and pains.

Operation principle:

Breo Temple and Eye Massager is used heat and pressure function help relax eye muscle and relieve minor eyes muscle aches and pains.

Product Structure

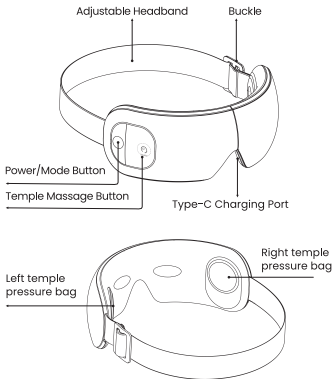



What's Included

1*Breo Temple and Eye Massager


1*USB charging Cable

1*User Manual

1*Storage Bag

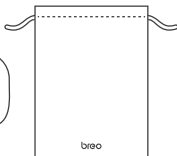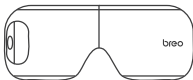

How to Use

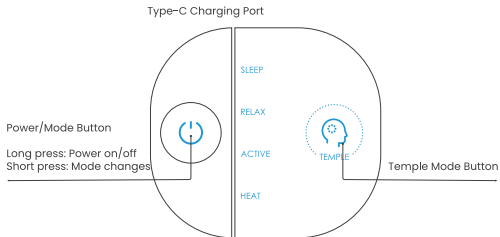


MODE	INDICATOR	MESSAGE
SLEEP	SLEEP light on	Soft massage with heat before sleep
RELAX	RELAX light on	For daily relaxation without heat
ACTIVE	ACTIVE light on	Stronger massage with heat
HEAT	HEAT light on	Heat Only
Left+Right Temple	SLEEP+RELAX+ACTIVE+ HEAT lights on	Left+Right Temple Massage
Left Temple	SLEEP+RELAX lights on	Left Temple Massage
Right Temple	ACTIVE+HEAT lights on	Right Temple Massage

How to Connect Breo App


1.Download the “breo” App

1.Search “breo” and download breo app from app store.

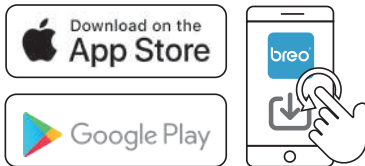


2.Scan here to connect with breo app.



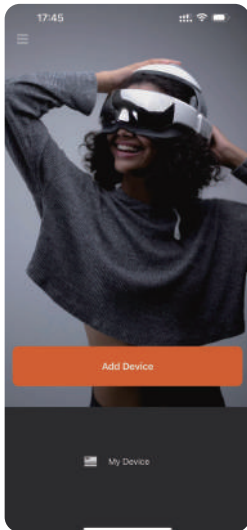

2.Add Device


- (1) Open the breo APP, log in App through your phone number. No need to log in repeatedly when entering the APP next time.
- (2) Long press "Power/Mode button" for 3 seconds to turn on the product. The massager must be turned on to connect to the APP.
- (3) Select "Add Device" after logging in, click product model "Breo_EM005_BLE" and connect the device.

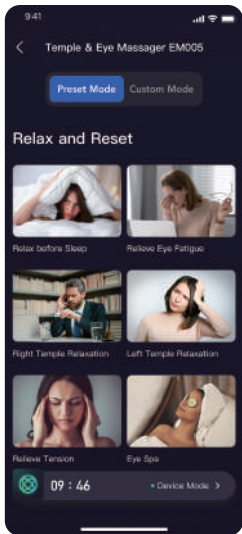

3.Preset Mode


Enter Breo Temple and Eye Massager user interface and tap on “Preset Mode”.

- (1) There are 6 preset modes for selection per your preference.
- (2) Select the mode you like and press “START” to start.
- (3) Press “PAUSE” to pause the mode you are using.

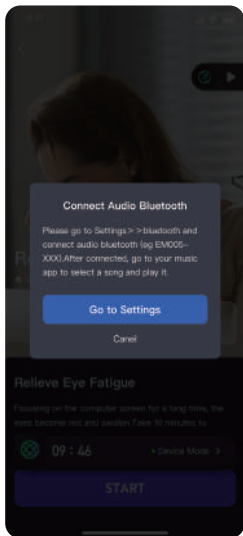

A prompt will appear on the interface indicating that audio Bluetooth needs to be connected.

Please go to Settings and connect audio Bluetooth "Breo_EM005"

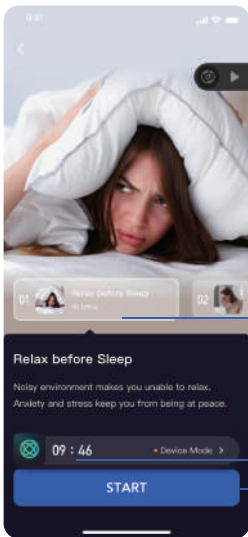

Green music icon:
Audio bluetooth connected successfully



Grey music icon:
Audio bluetooth connected unsuccessfully



Music playing,press it to pause music




Music pause,press it to start music

The current preset mode you are using.
You can swipe left or right to select a mode.

10 minute timer


Start/Pause the current preset mode

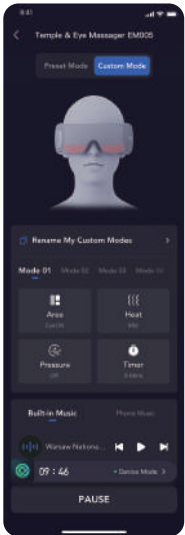

4.Custom Mode

Enter Breo Temple and Eye Messenger user interface and tap on "Custom Mode".

4.1 AREA

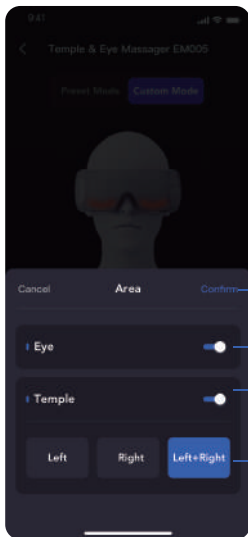


(4) After setting, press "Confirm" to save it.

(1) Tap on "Area" to set massage area.

(2) Eye or Temple area can be switched on/off individually.

(3) There are 3 options for the temple area: Left temple only-Right temple only-Left and right temple area.

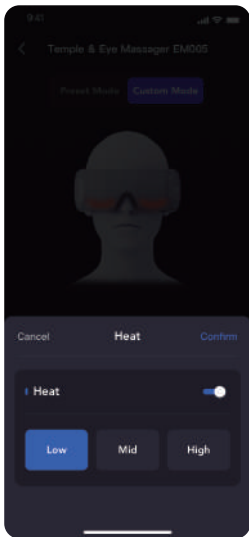


4.2 HEAT

(1) Tap on “Heat” to switch between low, mid and high level per your preference.


(2) Heat function can be switched on/off individually.

(3) After setting, press “Confirm” to save it.

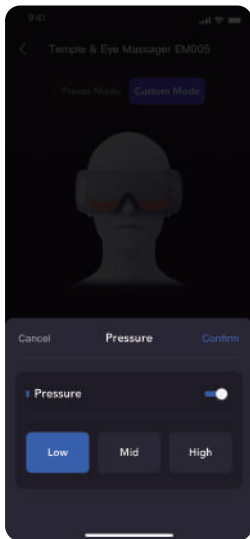

4.3 PRESSURE


- (1) Tap on "Pressure" to switch between low, mid and high level per your preference.
- (2) Pressure function can be switched on/off individually.
- (3) After setting, press "Confirm" to save it.

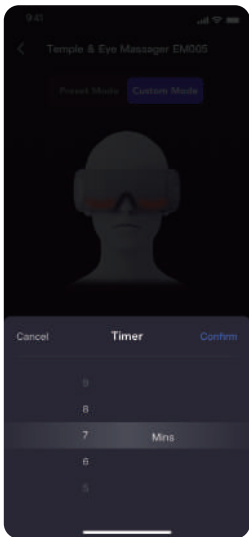

4.4 TIMER

(1) Tap on "Timer" to customize the working duration from 3mins to 20mins.


(2) After setting, press "Confirm" to save it.

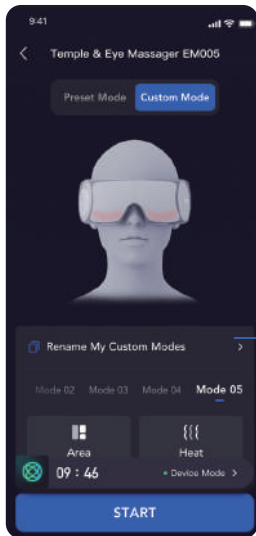

4.5 RENAME CUSTOM MODES

(1) You can save 5 custom modes and name them per your preferences.

(2) Press here to rename your custom modes.

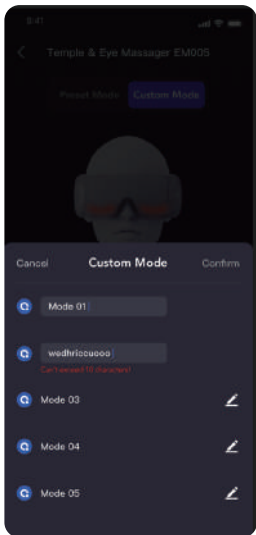

(3) You can rename the modes that you've set by tapping on the pen icon. Remember to press confirm to save your mode name.


(4) When you finish the settings, tap on the "Start/Pause" icon to start the massage.

5. Built-in Music / Phone Music


(1) Built-in Music : Click ► to play built-in music. Click ► to switch to the next built-in music song. Click ◀ to switch to the previous built-in music song. Click || to pause music.

(2) Phone Music:

- ① After choosing the Phone music , please go to your phone music app and select the music you like .
- ② Select the music you like, click ► to play it. And then go back to Breo App page.
- ③ If you want to switch to the next phone music song, please operate in your phone music app.

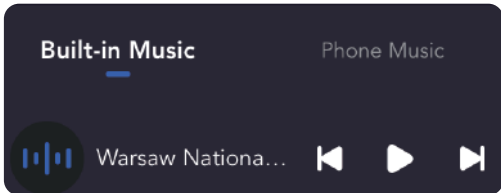

How to Wear

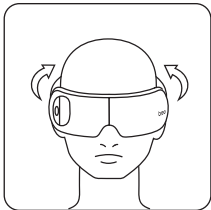

Align the cloth cover to your eyes and adjust the tightness of elastic band to the most comfortable position. Hold the power button for 3 seconds to turn on the eye massager.


Charging

1. Plug the USB charging cable into the socket of the charger.
2. Connect the charging cable into the charging port of the device.
3. Plug the charger into the power socket







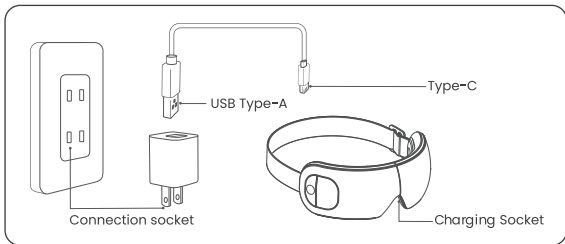


Operation	Indicator
Low Battery	SLEEP light flashes red
Automatic shutdown when battery is too low	SLEEP light flashes red three times and then shuts down.
Charing	SLEEP light flashes red
Charged full	SLEEP lights green
APP controlled	SLEEP+RELAX+ACTIVE+HEAT lights on

Caution:

1. It needs 2-3 hours for charging.
2. During the charging process, the product cannot be turned on.
3. If the product is charged while it is working, the product will automatically stop working.
4. When the product is not used for a long time, the battery will discharge slowly. In order to avoid battery damage due to low voltage for a long time, please charge the device for every three months.
5. Don't try to replace the rechargeable battery by yourself. Please contact the manufacturer for replacement. Self-disassembly and replacement the battery may cause damage to the main unit and battery.



Symbol description

	Keep Dry		Applied part of type BF
	Refer to instruction manual		Serial number
	Manufacturer information: The manufacturer name and address		Disposal in accordance with Directive 2002/96/EC (WEEE)
	Medical device		The device complies with part 15 of the FCC rules
	This device has not been tested for use in an MR environment and should not be used exposed to MR environments while patients are wearing the device. Keep it outside the MRI scanner room.		
RoHS	The device complies with compliance with RoHS Directive 2011/65/EU and its amendment directives 2015/863/EU		

Cleaning and Maintenance

Cleaning

1. Before cleaning, turn off the power and unplug the power adapter from the socket. Use a soft cloth with alcohol to clean the surface of the product.
2. Do not use corrosive cleaners such as benzene or diluents to clean the product, otherwise it will cause discoloration, surface aging, structural damage and other conditions.

Storage

1. If you do not use it for a long time, please cover it with a storage bag, and place the product in a dry and ventilated room to avoid sun exposure.
2. This product contains batteries, do not approach high temperature or put in the fire.

Troubleshooting

Problem	Possible causes	Take actions
Eye massager turn off automatically	1. The massage time was finished. 2. The Battery power is low.	1.You can long press Power button 3s to turn on the eye massager for another massage. 2.Please charge the eye massager.
The intensity is so strong	1. Elastic band is too tight. 2.May choose the hard intensity.	1.Please loosen the elastic band to the comfortable position. 2.Please choose the soft intensity on device or on APP.
The intensity is so weak	1. Elastic band is too loosen. 2.May choose the soft intensity.	1.Please tighten the elastic band to the comfortable position. 2.Please choose the hard intensity on device or on APP.
The device can not be charged	1.Please check power supply at home is normal. 2.The adapter may be overvoltage 5V, so it cannot charge the eye massager. 3.USB cable or adaptor may be damaged.	1.Please Please make sure that the power supply at home is normal . 2.Please use the adaptor which support 5V/1A to charge the eye massager. 3.Please use another adapter or USB cable to charge it.

Guarantee

This product comes with 1-year warranty from the date of purchase.

Warranty will be voided due to any of the following:

Damaged due to the improper usage or storage.

Damaged by attempting to repair or disassemble the device.

Period of warranty exceeded.

Product parameters

Product Name: Breo Temple and Eye Massager

Product Model: EM005

Charging Input: 5V=1A

Rated Power: 5W

Power Supply: 3.8V=2100mAh rechargeable lithium battery

Net Weight: 750g

Automatic power off: 10 minutes

Charging time: ≤3 hours

Working time: 70 minutes

Dimension: L220*W88*H112mm

FCC Regulatory Compliance

Contains FCC ID: PXU-BT-A1

This device complies with part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) this device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

FCC Warning Statement

Changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment. This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to Part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

Disposal

Notice:



At the end of the product life cycle, do not throw this product into the normal household garbage, but bring it to a collection point for the recycling of electronic equipment and battery.

Waste Electrical and Electronic Equipment may have potentially harmful effects on the environment.

Improper disposal may lead to the accumulation of harmful toxins in air, water and soil, which is harmful to human health. you have obligation to dispose of device correctly. Consult your municipal authority for information about disposal.

Customer Service

Please feel free to contact us if you have any questions or concerns.
Breco Customer Service Team guarantees a quick response and hassle-free solutions to any issue you may have within business hours.

Email: support@breocare.com

Web: us.breo.com

Tel: 1-866-546-3587

Working time :

PST 6:00am--11:00pm (Monday to Friday)

PST 7:00am--3:00pm (Saturday to Sunday)

Made in China



Shenzhen Breco Technology Co.,Ltd

19/F, Hisense Southern Building, No.1777 Chuangye Road, Nanshan District,
Shenzhen, Guangdong, CN, 518056

SIMPLY RELAX

