



© 2022 ALLSTAR INNOVATIONS.  
ALL RIGHTS RESERVED.  
CALMING HEAT IS A TRADEMARK OF ALLSTAR INNOVATIONS.  
MADE IN CHINA | PATENT PENDING  
ITEM# CWT 41003

CHBWON-3-IM1



SHARPER IMAGE®

## CALMING HEAT™ BACK WRAP

MODEL:CHBWON-3



NEED ASSISTANCE? PLEASE VISIT,  
[WWW.CALMINGHEATTIPS.COM](http://WWW.CALMINGHEATTIPS.COM)  
FOR ADDITIONAL INFORMATION ABOUT YOUR BACK WRAP.



**SHARPER IMAGE®**  
**CALMING HEAT™ BACK WRAP**



MODEL: CHBWON-3

## FCC CERTIFIED

---

This device complies with part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) This device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

Note: This product has been tested and found to comply with the limits for a Class B digital device, pursuant to part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This product generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this product does cause harmful interference to radio or television reception, which can be determined by turning the product off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the product and receiver.
- Connect the product into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

Changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

This Class B digital apparatus complies with Canadian ICES-003.

# FAQ

---

**How should I wash the Back Wrap?**

The Back Wrap is not machine washable and should only be spot cleaned with a damp cloth.

**Why doesn't my Back Wrap feel hot enough?**

Ensure that the product is preheated for 5 minutes with the heat setting on high.

**What should I do if my controller won't turn the Back Wrap on?**

Unplug the cord from the power source, wait 10 seconds, and plug the unit in again to reset it.

**What if I leave my Back Wrap on by accident?**

The Back Wrap comes with a 2-hour auto-off feature and a 15 minute auto-off for the vibration, and lumbar features.

# TABLE OF CONTENTS

---

<b>Important Safety Instructions</b>	<b>4</b>
<b>Reset and Instructions</b>	<b>6</b>
<b>Controller Features</b>	<b>7</b>
<b>Additional Waist Extension Feature</b>	<b>8</b>
<b>Cleaning Instructions</b>	<b>9</b>
<b>FAQ</b>	<b>10</b>
<b>FCC CERTIFIED</b>	<b>11</b>



# IMPORTANT SAFETY INSTRUCTIONS

## DO NOT DESTROY

### DANGER:

TO REDUCE RISKS OF BURNS, ELECTRIC SHOCK, FIRE, AND PERSONAL INJURY,  
THIS PRODUCT MUST BE USED IN ACCORDANCE WITH THE FOLLOWING INSTRUCTIONS:

1. READ ALL INSTRUCTIONS.
2. DO NOT USE WHILE SLEEPING.
3. DO NOT USE ON AN INFANT OR ON AN ANIMAL.
4. THIS PAD IS NOT TO BE USED ON OR BY AN INVALID, SLEEPING, HELPLESS OR UNCONSCIOUS PERSON, A PERSON WITH DIABETES, A PERSON WITH POOR BLOOD CIRCULATION OR INSENSITIVITY TO HEAT.
5. THIS APPLIANCE IS NOT TO BE USED BY CHILDREN. CLOSE SUPERVISION IS NECESSARY WHEN ANY APPLIANCE IS USED NEAR CHILDREN. CHILDREN SHOULD NOT PLAY WITH THE APPLIANCE.
6. DO NOT USE ON AREAS OF INSENSITIVE SKIN.
7. BURNS CAN OCCUR REGARDLESS OF CONTROL SETTING, CHECK SKIN UNDER PAD FREQUENTLY.
8. DO NOT USE THIS PRODUCT IF YOU HAVE SERIOUS CARDIOVASCULAR AND/OR RESPIRATORY DISEASE.
9. DO NOT USE THIS PRODUCT IF YOU HAVE A VERTEBRA, LUMBAR, SACRAL, VERTEBRAE, PELVIS, OR RIB FRACTURE.
10. DO NOT USE THIS PRODUCT IF YOU HAVE RHEUMATOID ARTHRITIS.
11. DO NOT USE THIS PRODUCT IF YOU HAVE ACUTE STAGE LUMBAR TRAUMA.
12. DO NOT USE THIS PRODUCT IF YOU HAVE PURULENT TUBERCULOUS SPONDYLITIS.
13. DO NOT USE THIS PRODUCT IF YOU HAVE GANGLION OR SKELETON LUMBAR.
14. DO NOT USE IN AN OXYGEN ATMOSPHERE.
15. DO NOT USE PINS OR OTHER METALLIC MEANS TO FASTEN THIS PAD IN PLACE.
16. DO NOT CRUSH PAD – AVOID SHARP FOLDS. PLACE PAD ON TOP OF AND NOT UNDER THE PART OF BODY NEEDING HEAT.
17. UNPLUG WHEN NOT IN USE. NEVER LEAVE APPLIANCE UNATTENDED, ESPECIALLY IF CHILDREN ARE PRESENT.
18. NEVER PULL THIS PAD BY THE SUPPLY CORD. DO NOT USE THE CORD AS A HANDLE. DO NOT PINCH CORD.

# CLEANING INSTRUCTIONS

### CAUTION:

DO NOT DRY CLEAN!  
DO NOT BLEACH!  
DO NOT USE WRINGER!

1. Disconnect pad from outlet.
2. Spot clean with damp rag or mild detergent
3. Lay flat to dry
4. Do NOT submerge in water
5. Do not plug in pad and pad connector until they are completely dry.

# PRODUCT SPECIFICATIONS


Pad Size:	52" x 9.6" (132.08 cm x 24.28 cm)
Power Rating:	Input: 100-240V 50/60Hz 0.35A Max Output: 5V 2000mA
Auto Shut-Off Timer	Heat 2 Hours Massage 15 Minutes
3 Massage Settings: Auto Off 15 Minutes	Low Speed Constant Med Speed Constant High Speed Constant
3 Heat Settings: Auto Off 2 Hours	Low Temp Range (~97F) Med Temp Range (~127F) High Temp Range (~135F)
6 Lumbar Settings: Auto Off 15 minutes	Low Inflate/Deflate Low Inflate and Hold Med Inflate/Deflate Med Inflate and Hold High Inflate/Deflate High Inflate and Hold

## ADDITIONAL WAIST EXTENSION FEATURE


This product includes a waist extension strap that will increase the overall length of the Back Wrap by 12 inches.

### To Install Extension Strap:

1. Attach the hook and loop portion of the extension strap to the inner facing portion of the Back Wrap (side facing towards body).

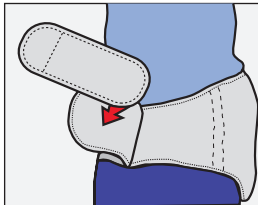

2. Place the partially assembled Back Wrap around your waist, and attach the extension strap to the hook and loop on the Back Wrap.

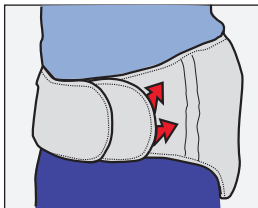


## AUTO-OFF FEATURE

This product includes an Auto-Off Feature: The heating function will turn off automatically after 2 hours of continuous use, and the vibration & traction feature will turn off after 15 minutes of continuous use.

**NOTE:** After a power outage, the Auto-Off feature automatically resets the timer.

19. DO NOT TUCK IN, TRAP, FOLD, CROSS, PINCH OR WRAP CONTROLLER OR CORD WITH HEATING PAD DURING USE.
20. NEVER LEAVE THIS PAD CONNECTED TO AN ELECTRICAL SUPPLY CIRCUIT WHEN NOT IN ACTUAL USE. WHEN NOT IN USE THIS PAD SHOULD BE KEPT IN A COOL, DARK, DRY PLACE – PREFERABLY IN THE ORIGINAL CARTON.
21. CAREFULLY EXAMINE PAD, COVER, CORD AND CONTROLLER BEFORE EACH USE. DISCARD THE PAD IF ANY SIGN OF DETERIORATION, WEAR OR DAMAGE. DO NOT TAMPER WITH THIS PAD IN ANY WAY. THERE ARE NO USER SERVICEABLE PARTS. IF PAD DOES NOT FUNCTION SATISFACTORILY, SEE WARRANTY OR CONTACT CUSTOMER SERVICE.
22. LOOP CORD LOOSELY WHEN STORING, TIGHT WRAPPING MAY DAMAGE CORD AND INTERNAL PARTS. ALLOW TO COOL BEFORE STORING AND DO NOT CREASE BY PLACING ITEMS ON TOP OF IT. STORE IN THE ORIGINAL PACKAGING, IN A COOL, DARK, DRY PLACE.
23. USE THIS PAD ONLY AT RATED VOLTAGE, AS MARKED ON THE PRODUCT.
24. DO NOT USE WITH GENERATORS, POWER CONVERTERS, OR INVERTERS.
25. DO NOT USE THIS PAD WITH LINIMENT, SALVE OR OINTMENT PREPARATIONS. SKIN BURNS COULD RESULT.
26. THIS APPLIANCE IS NOT INTENDED FOR MEDICAL USE IN HOSPITALS.
27. MODELS WITH A MASSAGE FUNCTION: DO NOT MASSAGE WHEN PAIN IS PRESENT, OR AREAS THAT ARE SWOLLEN, INFLAMED, COVERED WITH SKIN ERUPTIONS OR UNEXPLAINED PAIN. CONSULT A PHYSICIAN BEFORE USING ANY EQUIPMENT OF THIS TYPE.
28. USE ONLY WITH LISTED/CERTIFIED ITE POWER SUPPLY
29. DO NOT DRY CLEAN, DO NOT BLEACH, DO NOT IMMERSE IN WATER.
30. NEVER USE PAD IF THE PRODUCT FABRIC IS NOT INTACT.
31. READ AND FOLLOW ALL INSTRUCTIONS ON BOX OR PACKED WITH PAD BEFORE USING.
32. SAVE THESE INSTRUCTIONS.

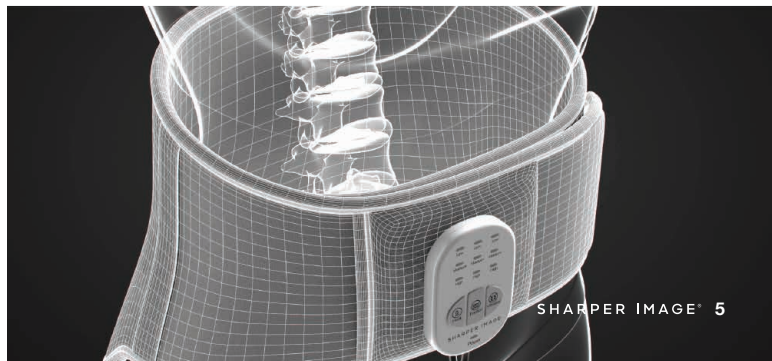

## INSTRUCTIONS

1. Plug the included USB adapter into any 120 volt AC electrical outlet. You can also power the Back Wrap with any 2.1A USB power source. Please note 1A USB power source will not achieve proper heat and vibration performance.
2. Press the Heat button 1X and allow product to preheat 5 minutes.  
Ensure the heat button is illuminated RED on high setting.
3. Wrap the heating pad around waist with the controller facing away from your body.
4. To turn either of the features off, cycle through the settings until LED light powers off.
5. Avoid crushing pad. Avoid sharp folds during use.



## RESET INSTRUCTIONS

If your heating pad unexpectedly shuts off (not including the Auto-Off setting feature), follow these steps:

1. Unplug the unit.
2. Wait 5 minutes.
3. Plug unit in and operate normally.

## CONTROLLER FEATURES

### Heat

Press 1x = High  
Press 2x = Med  
Press 3x = Low  
Press 4x = Off

### Vibration Massage


Press 1x = High Intensity Continuous  
Press 2x = Medium Intensity Continuous  
Press 3x = Low Intensity Continuous  
Press 4x = Off

**Note:** All 3 features can be turned on/off individually

### Lumbar

Press 1x = High Intensity:  
Inflate/Deflate  
Press 2x = High Intensity:  
Continuous Hold  
Press 3x = Medium Intensity:  
Inflate/Deflate  
Press 4x = Medium Intensity:  
Continuous Hold  
Press 5x = Low Intensity:  
Inflate/Deflate  
Press 6x = Low Intensity:  
Continuous Hold  
Press 7x = Off

**Note:** All 3 lumbar settings must deflate for 40 seconds. Please wait while the feature starts to run.

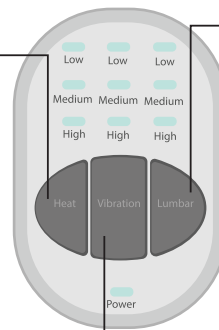

### LED Light Intensity Indicator

Low Medium High