


AquaSonic icon

Rechargeable Water Flosser

User Manual

Model: AS-ICON-FLOS

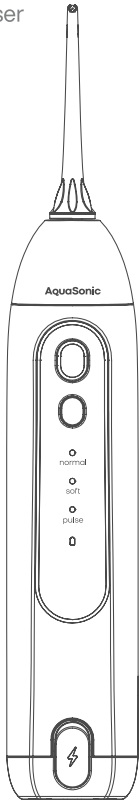

Before using this product, please read these instructions completely and save this manual for future use



IMPORTANT SAFEGUARDS

Before using this water flosser, please read the User Manual, and any other inserts that came with this product. The safety instructions provided here are intended to help you use this product correctly and safely, prevent you from injuring yourself and others, and avoid damage to property.

DANGER

- Keep out of reach of children. This device should not be used by toddlers, infants, or any people who are unable to operate it properly.
- Never alter this device, try to disassemble it, or repair it.
- Do not place this product where it can fall or be pulled into a bathtub, sink, or toilet. Do not submerge this product in water.
- If this product falls into water while charging, unplug it from the power outlet immediately.
- Do not use this product while bathing.

WARNINGS

- This product can be used by children age 8 and up, and by persons with reduced physical, sensory, or mental capabilities if they use it under supervision, or have been given adequate instruction by an adult who is responsible for their safety, and they understand the hazards involved.
- Children should not play with this device. Cleaning and user maintenance should not be done by children without supervision.
- This appliance is only to be used with the power supply unit provided with the appliance.
- The battery must be removed from this device before it is discarded. Ensure the product is unplugged and the battery power is fully drained before removing the battery. Dispose of the battery safely, according to local waste management regulations.
- Do not fill the tank with water warmer than 104°F / 40°C.
- If the charging cord is damaged, do not use it, to avoid a hazard.
- Keep the charging cord away from heated surfaces.
- Do not use this product without a nozzle.


1

- Do not use for any purpose other than oral cleaning.
- Keep nozzles away from small children. The nozzles are small enough to be swallowed and can become lodged in the throat.
- Certain types of mouthwash can damage this device by causing its housing or water reservoir to crack. Therefore, do not use mouthwash in the water reservoir. Use water only.
- People who experience tooth or gum pain may not be able to use this device and should consult their dentist.
- Do not press the nozzle too firmly against your teeth or gums. Doing so may injure your gums.

SAVE THESE INSTRUCTIONS

Questions? Contact support@aquasonic.com for assistance.

CONFIGURATION

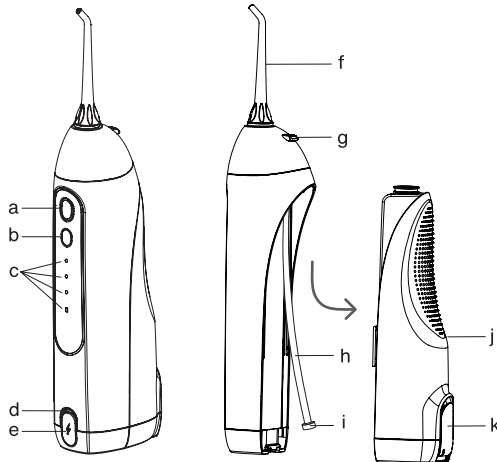

2

- a. On/Off button
- b. Mode button
- c. Mode/Charge indicators
- d. Charging port
- e. Charging port cover
- f. Nozzle
- g. Nozzle release button
- h. Nozzle suction tube
- i. Filter
- j. Water tank
- k. Water tank cap

Please note:

Upon initial unboxing, there may be some distilled water inside the water flosser that was used for its performance (water flow) quality assurance test. This is normal, and does not pose a health concern.

HOW TO USE

Charging the water flosser

1. Remove the charging port cover. Make sure the interior of the charging port is dry and the water flosser is turned off before charging.
2. Plug the charging cable into the charging port, and into a power outlet. (Fig. 1) The charging indicator light will blink, indicating that the device is charging.
3. When the battery is fully charged, the charging indicator light will turn off automatically. Unplug the charge cord, and re-insert the charging port cover.

Please note:

If the charging indicator light is blinking during or after the water flosser is used, that means the battery needs to be recharged. It takes about 4 hours to fully charge the battery.

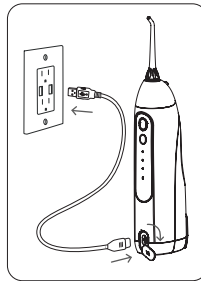


Fig. 1

3

Using the water flosser

1. Attach a nozzle to the water flosser. (Fig. 2)
2. Open the water tank cap, and fill it with cool or lukewarm water (Fig. 3), then close the water tank cap. Alternately, slide the water tank down off the unit (Fig. 4), fill it with water, then re-attach it to the unit. Water temperature must not exceed 104° / 40°C.
3. Press the Mode button to set your desired mode. You can also change the mode during flossing by pressing the Mode button.
 - a. **Normal** – standard flossing mode
 - b. **Soft** – gentler water pressure for sensitive teeth and gums
 - c. **Pulse** – pulsed water stream for massaging gums
4. Place the nozzle inside the mouth in an upright position. Gently close the mouth around the nozzle, and lean over the sink so water can run out. (Fig. 5)
5. Press the ON/OFF button to start flossing. Direct the water jet stream at a 90-degree angle to the teeth and gums. Move the stream slowly across the teeth. Clean between teeth and around any braces, orthodontia, or dental work. (Fig. 6)
6. Once water flossing is complete, press the ON/OFF button to turn off the device.
7. To change the nozzle, press the nozzle release button, and pull the nozzle up. (Fig. 7)

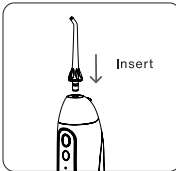


Fig. 2

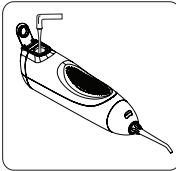


Fig. 3

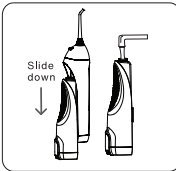


Fig. 4

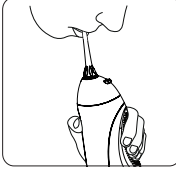


Fig. 5

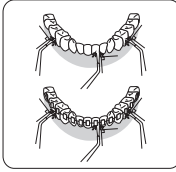


Fig. 6

4

After using the water flosser

1. Open the water tank cap and pour out any water that remains in the tank.
2. Press the ON/OFF button to turn on the device and drain any water from within the mechanism.
3. Press ON/OFF again to turn the device off, then close the water tank cap.
4. Wipe the device with a dry cloth. (Fig. 8)

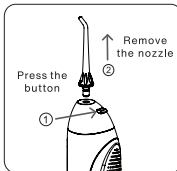


Fig. 7

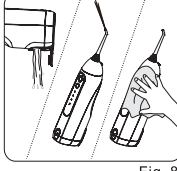


Fig. 8

Please note:

To prevent bacteria from propagating, make sure all water is drained from within the unit.

MAINTENANCE

Use only water or a neutral cleaning agent to clean this product. Never use caustic or abrasive agents, such as vinegar or lime scale removers, as they may damage the device. Do not wash using water that's over 122°F / 50°C. Keep the device away from heat, and avoid direct sunlight.

Main unit

1. To keep the device clean, wipe off any debris with a soft, damp cloth.
2. Do not submerge the main unit in water to clean it.
3. Dry the unit after each use and after cleaning.

Water tank cap

Wipe off any grime with a cloth.

Nozzle

1. Wash with water and dry with a soft cloth.
2. Replace the nozzle every 6 months.

5

Water suction tube

1. Wash with water and dry with a soft cloth.
2. Do not bend, pull or twist the water suction tube.

Tank

1. Wash with water and dry with a soft cloth.
2. Dry the tank when storing for a long period of time.

TROUBLESHOOTING

Problem	Possible cause	Solution
The unit doesn't work	You have just purchased the product or it has not been used for more than 3 months	Charge the main unit at least 4 hours
	There is dirt stuck in the nozzle	Replace the nozzle
	The nozzle is deformed	Replace the nozzle
Water pressure is very low	There is dirt stuck in the filter	Clean the filter
	The battery has reached the end of its life cycle	Contact service center for support
	Charging time was insufficient	Charge the main unit at least 4 hours
The unit operates only for a short period of time	The tank is empty	Fill up the tank
	The appliance is tilted too much	Use with the appliance in a vertical position
	You've just purchased the product or it has been left unused for a period of time	Pull out the nozzle then flush nozzle hole with water
Water does not come out		

6

DISPOSAL

At the end of this product's life cycle, remove the rechargeable battery from within the handle by following the instructions below. Protect your eyes, hands, fingers and your work surface while removing the battery, then dispose of it according to local waste management requirements.

Removing and recycling the battery

1. Slide the water tank off the unit.
2. Remove the rubber feet from the bottom of the handle, then use a screwdriver to release the bottom panel. (Fig. 9)
3. Remove the top cover using a screwdriver and unscrew the 2 screws to release the inner part. (Fig. 10)
4. Pull out the interior mechanism, and unscrew the screws. (Fig. 11)
5. Cut the battery line and remove the batteries. (Fig. 12) If you have disposal or recycling questions, please contact your local waste management office.

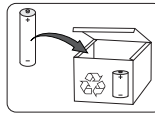


Fig. 9

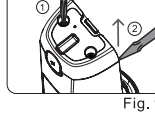


Fig. 10

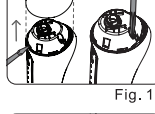


Fig. 11

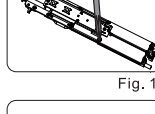


Fig. 12

DISCLAIMER

Users of this device do so at their own risk. Neither AquaSonic nor its retailers assume any responsibility or liability for any injury or damages, physical or otherwise, resulting, directly or indirectly, from the use of this device. Further, AquaSonic reserves the right to revise this publication and to make changes from time to time in the contents thereof without obligation to notify any person of such revision or changes. Model may be changed for improvements without notice.



Printed in China

7