

折叠尺寸: 58x190mm (不含3mm出血) 展开尺寸: 238x196 (含3mm出血) 纸张: 120g双胶纸 颜色: C0 M0 Y0 K70

## AquaSonic icon

### USER MANUAL



Model: AS-ICON

Before using this product, please read these instructions completely and save this manual for future use.

#### IMPORTANT SAFEGUARDS

Before using this toothbrush, please read the User Manual, and any other inserts that came with this product. The safety instructions provided here are intended to help you use this product correctly and safely, and to prevent you from injuring yourself or others.

#### DANGER-VOLTAGE

- Voltage converters do not guarantee compatibility.
- Do not force the plug into an outlet.
- If the plug does not fit easily into the outlet, discontinue use.
- Always unplug this product immediately after using.

#### MEDICAL WARNINGS


- If you have had oral surgery or gum surgery within the past two months, consult your dentist before using this product.
- Consult your dentist if excessive bleeding occurs while using this product or if minor bleeding continues to occur after one week of use.
- Consult your physician prior to using this toothbrush if you have any medical concerns.
- Consult your physician or implant device manufacturer prior to using this product if you have questions concerning the combined use of this product with a pacemaker or other implanted medical device(s).

#### SAFETY WARNINGS

- This toothbrush should not be used on children, or persons with reduced physical, sensory, or mental capacities unless supervised by a responsible adult party. Cleaning and user maintenance should not be done by children without supervision.
- Children should be supervised to ensure they do not play with this toothbrush.
- This toothbrush should only be used for its intended purpose as described in this user manual. It is designed to clean teeth and gums only. Do not use any attachments or brush heads other than those recommended by the manufacturer.
- This toothbrush is a personal care device and is not intended to be shared or used by more than one person.
- If your toothpaste includes peroxide, baking soda, or other bicarbonates which are common in whitening toothpastes, please ensure that you rinse your brush head and handle thoroughly between uses, to avoid cracking the plastic.
- Do not use this toothbrush if it is not functioning properly. If it has been dropped, damaged, or dropped into water, discontinue use immediately.
- Keep this product away from heat and open flame.
- If you experience any discomfort or pain in the mouth as a result of using this toothbrush, discontinue use and contact your dentist.
- Do not attempt to remove the battery.
- Do not clean this toothbrush in the dishwasher.
- Never drop or insert any object into any opening in the toothbrush.

**SAVE THESE INSTRUCTIONS.**

#### CONFIGURATION

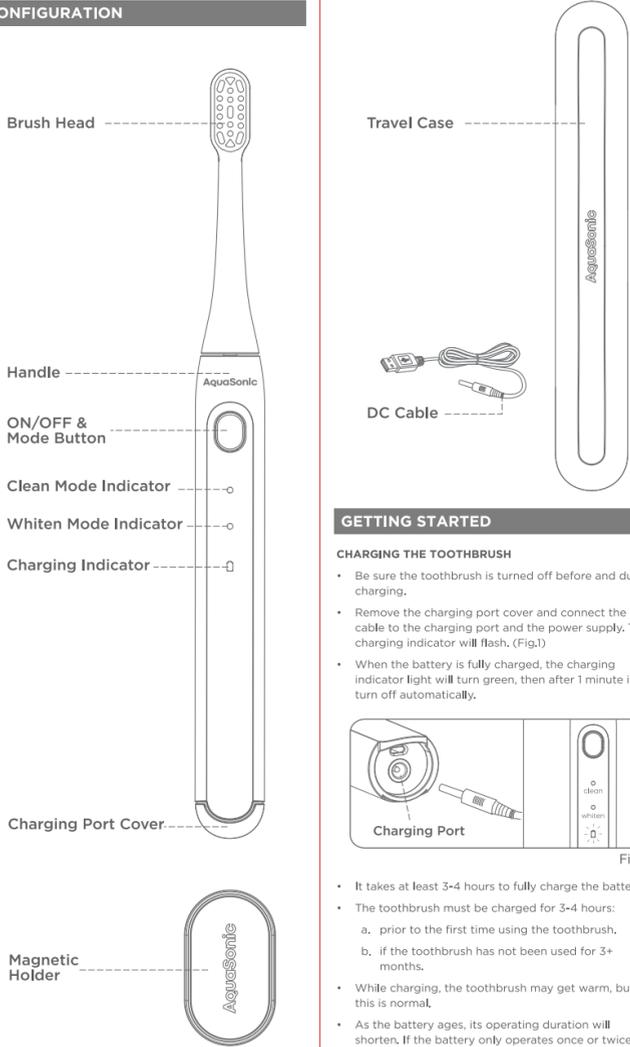

Brush Head

Travel Case

Handle

ON/OFF & Mode Button

Clean Mode Indicator

Whiten Mode Indicator

Charging Indicator

Charging Port Cover


Magnetic Holder

DC Cable

#### GETTING STARTED

##### CHARGING THE TOOTHBRUSH

- Be sure the toothbrush is turned off before and during charging.
- Remove the charging port cover and connect the DC cable to the charging port and the power supply. The charging indicator will flash. (Fig.1)
- When the battery is fully charged, the charging indicator light will turn green, then after 1 minute it will turn off automatically.

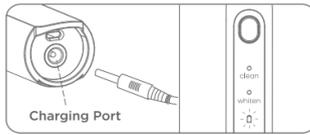

Charging Port

Fig.1

- It takes at least 3-4 hours to fully charge the battery.
- The toothbrush must be charged for 3-4 hours:
  - prior to the first time using the toothbrush.
  - if the toothbrush has not been used for 3+ months.
- While charging, the toothbrush may get warm, but this is normal.
- As the battery ages, its operating duration will shorten. If the battery only operates once or twice after being fully charged, it is nearing the end of its life cycle. Contact support@aquasonic.com for support.

#### HOW TO USE

- Slide the brush head onto the handle and push down firmly. (Fig.2) There should be a small gap between the base of the brush head and the handle.
- Wet the bristles and apply toothpaste.
- Place the brush head in your mouth.
- Lightly press the button on the handle (Fig.3) and start brushing your teeth. A mode indicator will light up:
  - Clean - standard mode, for an everyday clean.
  - Whiten - polishing mode, to help reduce the look of surface stains.
- Brush evenly throughout the mouth, in 4 sections (Fig.4): 1-outside upper, 2-outside lower, 3-inside upper, 4-inside lower.
- Holding the bristles at a 45-degree angle to the gumline, guide the brush head slowly from tooth to tooth, pausing for a few seconds on each tooth. Brush the gums and the chewing surfaces, as well. Do not press too hard or scrub your teeth - simply let the toothbrush do all the work.
- During brushing, the vibration will momentarily pause every 30 seconds, indicating to move to the next quadrant of the mouth. After 2 minutes, the toothbrush will shut off automatically.
- If you press the button while the toothbrush is on, it will change modes, or turn off. (Fig. 5 & Fig.6)

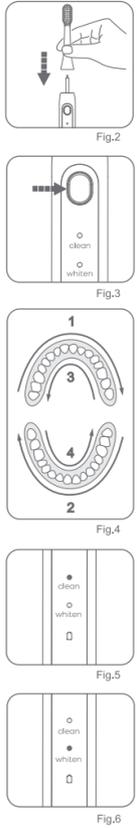


Fig.2

Fig.3

Fig.4

Fig.5

Fig.6

NOTE:

- 2 minutes is generally long enough to clean your teeth and gums properly. If you'd like, you may brush for longer than 2 minutes, for an even more thorough clean.
- It is common to experience slight bleeding when switching from a manual toothbrush to an electric toothbrush for the first time. This should stop after a few days, when your gums get used to the toothbrush. If bleeding persists, please consult your dentist.

#### CLEANING & MAINTENANCE

- Turn on the toothbrush and rinse the brush head under a tap. (Fig.7)
- Turn off the toothbrush, detach the brush head, and rinse the inside of the brush head with cold or lukewarm water. (Fig.8)
- Turn on the toothbrush handle and rinse the drive shaft under a tap to remove toothpaste residue. (Fig.9)
- Put the brush head back on the handle.
- Wipe the entire toothbrush with a soft damp cloth or paper towel.

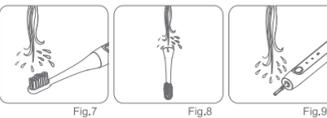


Fig.7

Fig.8

Fig.9

NOTE:

For optimal results, replace the brush head every 3 months, or sooner if the bristles become worn and start to bend outwards.

#### TROUBLESHOOTING

PROBLEM	POSSIBLE CAUSE	SOLUTION
The toothbrush doesn't work.	You've just purchased the product or it has been left unused for more than 3 months.	Charge the toothbrush handle for at least 3-4 hours.
Mode indicator lights up, but the toothbrush doesn't work.	Battery is low.	Charge the toothbrush handle for at least 3-4 hours.
Weak or diminished oscillations.	You're pressing the toothbrush up against your teeth too hard.	Place the bristles gently against your teeth.
The toothbrush only works for a few minutes, even after charging.	The battery has reached the end of its life cycle.	Contact AquaSonic for support.
	Charging time was insufficient.	Charge the toothbrush handle for at least 3-4 hours.

#### DISPOSAL

AquaSonic Icon is powered by a lithium-ion battery, which can be recycled. At the end of this product's life, always remove the battery before discarding the handle, and recycle or dispose of the battery according to local waste management requirements. If you have any questions, please contact your local waste management service provider.

#### REMOVING AND RECYCLING THE BATTERY

- Before removing the battery, make sure its charge is completely drained.
- With pliers, clamp the handle approximately 1-1.5 inches above the base of the handle. Use a flathead screwdriver to push the post of the bracket until it is released. (Fig.10)
- Use an electric iron to melt the tin that connects the battery and PCB. Take out the battery (Fig.11) and recycle it (Fig.12)

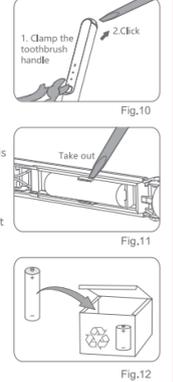


Fig.10

Fig.11

Fig.12

#### DISCLAIMER

Users of this device do so at their own risk. Neither AquaSonic nor its retailers assume any responsibility or liability for any injuries or damages, physical or otherwise, resulting, directly or indirectly, from the use of this device. Further, AquaSonic reserves the right to revise this publication and to make changes from time to time in the contents thereof without obligation to notify any person of such revision or changes. Model may be changed for improvements without notice.

