

Brookstone®

TheraSqueeze™ Foot and Calf Massager

BRING TIRED CALVES & FEET BACK TO LIFE

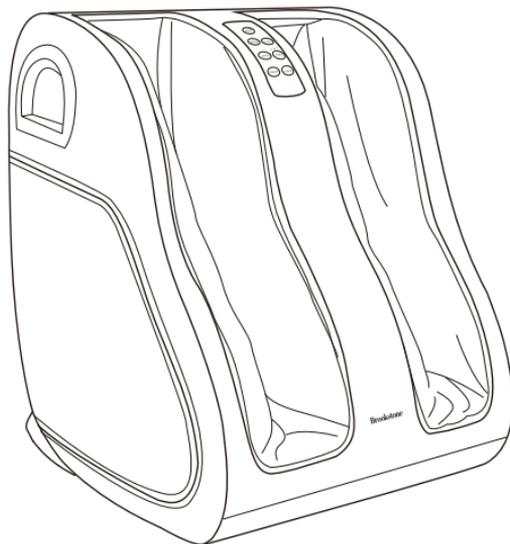


TABLE OF CONTENTS

Warnings and Caution	2
Location of Parts	7
Location of Controls	8
Preparation Before Using Power On/Off	9
Using the Massager	10
Product Care and Storage	11
Specifications	11
Troubleshooting Guide.....	12
Warranty.....	13

CAUTION

- TO REDUCE THE RISK OF ELECTRIC SHOCK, DO NOT REMOVE COVER. THERE ARE NO SERVICEABLE PARTS INSIDE.
- TO REDUCE THE RISK OF FIRE OR ELECTRIC SHOCK, DO NOT EXPOSE THE UNIT TO RAIN OR MOISTURE.



The lightning flash with arrow-head symbol within an equilateral triangle is intended to alert the user to the presence of uninsulated “dangerous voltage” within the unit’s enclosure that may be of sufficient magnitude to constitute a risk of electric shock.



The exclamation point within an equilateral triangle is intended to alert the user to the presence of important operating and maintenance (servicing) instructions in the literature accompanying the unit.

IMPORTANT SAFETY INSTRUCTIONS

All of the safety and operating instructions should be read, adhered to and followed before the unit is operated.

SAVE THESE INSTRUCTIONS

When using an electrical appliance, basic precautions should always be followed. Read all instructions before using this unit.

DANGER

To reduce the risk of electric shock, burns, fire or injury:

1. Always unplug this massager from the electrical outlet immediately after use and before cleaning.
2. Do not use while bathing or in a shower.
3. Do not place or store unit where it can fall or be pulled into a tub or sink.
4. Do not place, drop or submerge in water or other liquid.
5. Do not reach for unit that has fallen into water. Unplug it immediately.
6. Care should be taken so that objects do not fall and liquids are not spilled onto the unit.
7. Never use pins or other metallic fasteners with this appliance.
8. Carefully examine the covering before each use. Discard the appliance if there are any signs of deterioration such as cracking or blistering.
9. Keep dry—do not operate in a wet or moist environment.
10. Do not use wet hands to connect or disconnect the power plug from the electrical outlet.

WARNING - To reduce the risk of burns, fire, electric shock, or injury to persons: Read all instructions before using

1. The unit should never be left unattended when plugged in. Unplug from outlet when not in use and before putting on or taking off parts.
2. Do not operate under blanket or pillow. Excessive heating can occur and cause fire, electric shock, or injury to persons.
3. Close supervision is necessary when this unit is used by, on, or near children, invalids, or disabled persons.

4. Use this unit only for its intended use as described in this instruction manual. Do not use attachments not recommended by Brookstone.
5. Never operate this unit if it has a damaged cord or plug, if it is not working properly, if it has been dropped or damaged, or if it has been dropped into water. Contact Brookstone Customer Service at 866-336-3199.
6. Consult your physician before using this unit.
7. The unit is not intended for self-treatment of conditions that should be managed by a qualified health care provide.
8. The unit is designed as a non-professional massager to soothe tired and aching muscles. It is not intended to be a substitute for medical attention. If your ailments and symptoms persist or worsen after use, contact your physician immediately.
9. **IMPORTANT:** Any individual who may: be pregnant or recently gave birth, has a heart condition or a pacemaker, suffers from diabetes or other sensory impairment, phlebitis and/or thrombosis, joint dysfunction, inflammation or injuries, weak bones, osteoporosis, a high fever, an abnormal or curved spine, is at an increased risk of developing blood clots, has pins/screws/artificial joints or other medical devices implanted in his/her body, is being treated for a medical condition or physical injury or is under a doctor's care should consult with a medical professional before using a massaging device designed for home use.
10. Do not carry this unit by power cord or use power cord as a handle.
11. Keep the cord away from heated surfaces.
12. Never operate the unit with the air openings blocked. Keep the air openings free of lint, hair and the like.
13. Never drop or insert any object into any opening.
14. Do not use outdoors.
15. Do not operate where aerosol spray products are being used or where oxygen is

being administered.

16. To disconnect, turn all controls to the off position, then remove plug from outlet.
17. Unplug this unit during lightning storms or when unused for long periods of time.
18. Protect the power cord from being walked on or pinched, particularly at plug outlets, convenience receptacles and the point where it exits the unit.
19. Do not allow cord to touch hot surfaces. Wrap power cord loosely around the unit when storing.
20. Do not Crush - Avoid sharp folds.
21. Do not operate for more than 20 minutes continuously.
22. Do not perform continuous massage on the same spot of your body for more than 5 minutes at a time or this may result in excessive stimulation and may have adverse effects.
23. The unit should be situated away from direct sunlight or heat sources such as radiators, electric heaters, heat registers, stoves or other units (including amplifiers) that produce heat. Avoid placing on top of stereo equipment that radiates heat.
24. Never block the air openings of the unit with materials such as clothing, plastic bags or papers, or place it on a soft surface such as a bed or couch, where the air openings may be blocked.
25. Do not overload the electrical outlet. Use your power source as instructed.
26. To avoid the risk of electric shock, do not disassemble or attempt to repair the unit. Incorrect repair can cause risk of electric shock or injury to persons when unit is used.
27. Do not operate in the presence of explosive and/or flammable fumes.
28. Never remove the plug from the socket by pulling the power cord.
29. Do not use in an excessively humid environment, such as the bathroom.

30. Do not use on a wet body. This may result in a short circuit or cause electric shock.
31. Do not use if the fabric covers are torn or if parts are missing.
32. Do not exert pressure when using the unit. This may result in injury.
33. Turn the unit off and unplug it immediately if:
 - liquid is spilled on it.
 - you feel severe abnormal pain.
 - you detect any malfunction or any other abnormal condition.
 - there is a power failure or lightning.
34. Do not use the unit if the electrical outlet is excessively loose.
35. You may experience itching on your feet during and after the massage. This is normal and is a result of improved blood circulation.
36. Do not use massager in close proximity to loose clothing or jewelry.
37. Do not stand on the massager. Use only while seated.
38. This is designed for household use only and not commercial purposes.
39. Read instructions before using.
40. Never remove cover.

IMPORTANT

This product has a polarized plug (one blade is wider than the other).

As a safety feature, this plug will fit in a polarized outlet only one way.


If the plug does not fit fully in the outlet, reverse the plug. If it still does not fit, contact a qualified electrician. Do not attempt to bypass this safety feature.

SAVE THESE INSTRUCTIONS

LOCATION OF PARTS


Massage Unit

1. Control Panel
2. Roller Massage for Foot and Calf
3. Handles
4. Dynamic Squeeze™ Massage for Foot and Calf
5. Fabric Cover
6. Kickstand (not shown)
7. Power Cord
8. Power Switch (not shown)

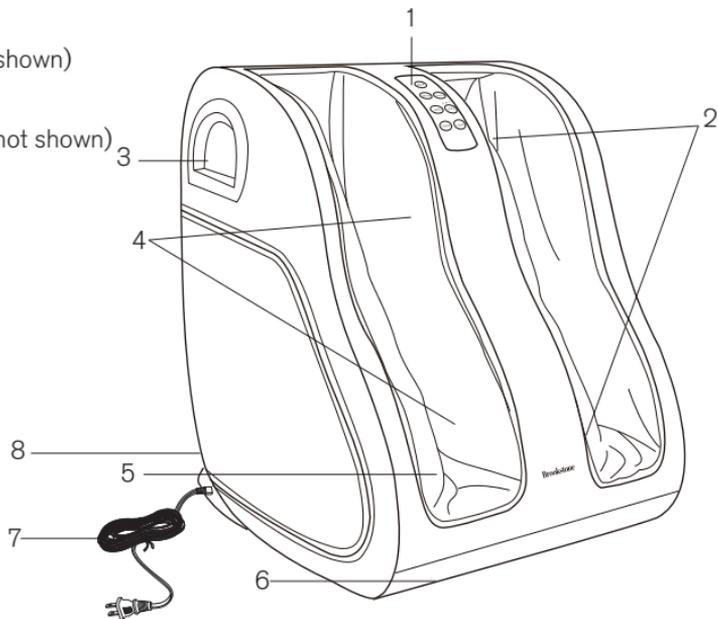

LOCATION OF CONTROLS

1. AUTO 1
 - Automatic program 1
2. AUTO 2
 - Automatic program 2
3. CALF
 - Calf automatic program
4. FOOT
 - Foot automatic program
5. INTENSITY (Low, Medium, High)
 - Adjustable for Dynamic Squeeze™ massage
6. ROLLER
 - Turn foot roller on/off
7. POWER BUTTON
 - Power on/off

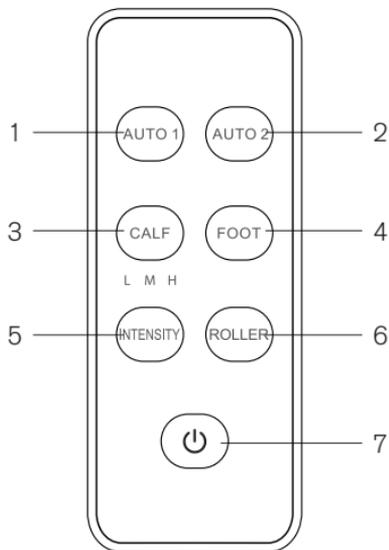

PREPARATION BEFORE USING POWER ON/OFF

1. Plug the power cord into a 110V power source.
2. Turn the power switch (located at the rear corner on the back of the unit) to ON.
3. Press power button on the control panel to turn the unit on.
4. Press power button again to turn the unit off.

USING THE MASSAGER

1. Plug the unit into a properly grounded wall outlet.
2. Press the MAIN POWER SWITCH on the back of the unit to on.
3. Press the POWER BUTTON to turn on the unit.
4. Press the AUTO BUTTON to turn the automatic program for Foot and Calf on/off.
5. Press the FOOT BUTTON to turn the automatic program for Foot on/off.
6. Press the CALF BUTTON to turn the automatic program for Calf on/off.
7. Press the INTENSITY BUTTON to switch between LOW, MEDIUM, HIGH air pressure.
8. Press the ROLLER BUTTON to turn Foot Rollers on/off.
9. The AUTO, FOOT, CALF programs each include different roller action cycles with air pressure massage. A 15-minute massage is recommended.
The unit will automatically turn off after 20 minutes of continuous operation.

Note: The massage may tickle, which is normal and a sign of blood circulation improvement.

10. To turn the unit off, press the POWER BUTTON.
11. If the massager will not be used for an extended period of time, press the MAIN POWER SWITCH to OFF.

IMPORTANT: Do not stand with your feet in the foot beds. Do not sit on the unit or press hard against the foot beds, as this may cause the motor to labor and overheat.

PRODUCT CARE AND STORAGE

- Turn off and unplug the unit before cleaning.
- Do not use volatile liquids such as benzene, thinner or gasoline for cleaning.
- Clean the main unit with a damp cloth or a neutral cleanser and then wipe with a dry cloth.
- Unzip the top layer of fabric cover from the unit and machine wash. Tumble dry on low.
IMPORTANT: keep the inner fabric cover dry before reattaching top layer to the main unit.
- Do not machine wash the inner fabric cover.
- If you don't use the unit for an extended time, roll up the power cord and keep the unit away from environments where it's humid and dusty.
- Do not expose the unit to heaters and direct sunlight.

SPECIFICATIONS

Voltage	AC 110/60Hz
Power Consumption.....	80W
Automatic Timer	20 minutes
Net Weight	27.6 lbs
Dimensions	17.5"l x 14.5"w x 19.5"h

TROUBLESHOOTING GUIDE

ISSUE

POSSIBLE CAUSE/RESOLUTION

The unit stops during operation.

The 20-minute auto shutoff is engaged. Restart massager using power button.

The thermal resetting protection device is engaged. This is an added safety feature to protect the unit from overheating. If the unit stops and fails to restart as a result of overheating, turn the main power off and do not operate the unit for at least 30 minutes. Failure to turn off the unit may cause it to resume operation unexpectedly when the device is cooled.

The main power switch or the power button on the control panel was pressed accidentally.

The unit does not start.

Ensure the power cord is connected to a properly working electrical outlet.

ONE (1) YEAR LIMITED WARRANTY

Brookstone® warrants this product against defects in materials and/or workmanship under normal use for a period of ONE (1) YEAR from the date of purchase by the original purchaser (“Warranty Period”). If a defect arises and a valid claim is received within the Warranty Period, at its option, Brookstone will either 1) repair the defect at no charge, using new or refurbished replacement parts, or 2) replace the product with a new product that is at least functionally equivalent to the original product, or 3) provide a store credit in the amount of the purchase price of the original product. A replacement product or part, including a user-installable part installed in accordance with instructions provided by Brookstone, assumes the remaining warranty of the original product. When a product or part is exchanged, any replacement item becomes your property and the replaced item becomes Brookstone’s property. When a store credit is given, the original product must be returned to Brookstone and becomes Brookstone’s property.

Obtaining Service: To obtain warranty service, call Brookstone Limited Warranty Service at 1-800-292-9819. Please be prepared to describe the product that needs service and the nature of the problem. A purchase receipt is required. All repairs and replacements must be authorized in advance. Service options, parts availability and response times will vary. You are responsible for delivery and the cost of delivery of the product or any parts to the authorized service center for replacement, per our instructions.

Limits and Exclusions: Coverage under this Limited Warranty is limited to the United States of America, including the District of Columbia and the U.S. Territories of Guam, Puerto Rico, and the U.S. Virgin Islands. This Limited Warranty applies only to products manufactured for Brookstone that can be identified by the “Brookstone” trademark, trade name, or logo affixed to them or their packaging. The Limited Warranty does not apply to any non-Brookstone products. Manufacturers or suppliers other than Brookstone may provide their own warranties to the purchaser,

but Brookstone, in so far as permitted by law, provides these products “as is.” This warranty does not apply to: a) damage caused by failure to follow instructions relating to product’s use or the installation of components; b) damage caused by accident, abuse, misuse, fire, floods, earthquake or other external causes; c) damage caused by service performed by anyone who is not a representative of Brookstone; d) accessories used in conjunction with a covered product; e) a product or part that has been modified to alter functionality or capability; f) items intended to be periodically replaced by the purchaser during the normal life of the product including, without limitation, batteries or light bulbs; g) any product sold “as is” including, without limitation, floor demonstration models and refurbished items; or h) a product that is used commercially or for a commercial purpose.

BROOKSTONE SHALL NOT BE LIABLE FOR INCIDENTAL OR CONSEQUENTIAL DAMAGES RESULTING FROM THE USE OF THIS PRODUCT, OR ARISING OUT OF ANY BREACH OF THIS WARRANTY. TO THE EXTENT PERMITTED BY APPLICABLE LAW, BROOKSTONE DISCLAIMS ANY AND ALL STATUTORY OR IMPLIED WARRANTIES, INCLUDING, WITHOUT LIMITATION, WARRANTIES OF MERCHANTABILITY, FITNESS FOR A PARTICULAR PURPOSE AND WARRANTIES AGAINST HIDDEN OR LATENT DEFECTS. IF BROOKSTONE CANNOT LAWFULLY DISCLAIM STATUTORY OR IMPLIED WARRANTIES, THEN TO THE EXTENT PERMITTED BY LAW, ALL SUCH WARRANTIES SHALL BE LIMITED IN DURATION TO THE DURATION OF THIS EXPRESS WARRANTY.

Some states disallow the exclusion or limitation of incidental or consequential damages or how long an implied warranty lasts, so the above exclusions or limitations may not apply to you. This warranty gives you specific legal rights and you may also have other rights, which vary from state to state.

Find thousands more
great ideas online

Brookstone®

314296

© Brookstone 2016

Merrimack, New Hampshire, USA 03054 ▪ 800-846-3000 ▪ Brookstone.com