



### Good Brushing Technique:

Proper brushing takes at least two minutes, but most adults don't even come close to brushing that long. To make it easier, the Dazzlepro Elite and Advanced Sonic Toothbrush models have been equipped with a Two Minute Smart Auto Timer and Area Change Alert which reminds you to brush each quadrant of your mouth for 30 seconds and then automatically shuts off once you have achieved the proper brushing time.

To experience the superior cleaning power of your Dazzlepro toothbrush, simply place the bristles gently against your teeth and let the sonic power go to work. Do not press too hard or scrub, simply let the brush do the work for you. Pay extra attention to the gum-line, hard-to-reach back teeth and areas around crowns, fillings or other dental work.

For each 30 second interval concentrate on thoroughly cleaning each of the following areas of your mouth:

- Clean the outer surfaces of your upper teeth, then outside your lower teeth
- Clean the inner surfaces of your upper teeth, then inside your lower teeth
- Clean all the chewing surfaces
- For fresher breath, be sure to brush your tongue, too

### How Important is the Toothpaste I Use?

It is important that you choose a toothpaste that's right for your oral needs. A wide selection of toothpastes are available and designed for nearly every condition, so ask your dental professional which toothpaste is right for you.



## TWO-YEAR WARRANTY

This Dazzlepro product is warranted to be free of manufacturer's defects in materials or workmanship for two years from the date of purchase. Damage or wear resulting from an accident, misuse, abuse, commercial use, or unauthorized adjustment and/or repair are not covered under this warranty.

Should this product require warranty service (or replacement at our discretion) please contact client service to obtain a Return Merchandise Authorization number (RMA) and return instructions, proof of purchase is required. Products returned without a Dazzlepro generated RMA number will not be accepted and the sender will not receive a refund, replacement, or repaired product.



info@dazzlepro.com  
1.888.358.0330

There are no express warranties except as listed above. This warranty gives you specific legal rights and you may have other rights which vary from state to state.

DO NOT RETURN TO RETAILER, PLEASE CONTACT US DIRECTLY FOR WARRANTY SERVICE OR REPAIR, THANK YOU.

Thank you for choosing  
**Dazzlepro**  
in the spirit of welcoming you  
we would like to offer ...

**50%  
OFF**

Your next purchase!  
Just enter this exclusive code

**UVCASE50**

at checkout when you visit  
dazzlepro.com

### Safety & Precautions

Always follow basic safety precautions when using electrical devices. The Dazzlepro Elite Sonic Toothbrush has a waterproof handle, it is electrically safe and specifically designed for use in the bathroom. If the device fails to function, do not attempt to dismantle, exchange or manipulate any part of it to avoid the risk of electric shock. This product is not intended for use by children under the age of three.

© 2015, Dazzlepro. All Rights Reserved

**Dazzlepro**<sup>TM</sup>  
UV SANITIZING CHARGING CASE



# Welcome to Dazzlepro™

Please take a moment to read these instructions thoroughly so that you may enjoy all the features that make the UV Sanitizing Charging Case the perfect accessory for your Dazzlepro toothbrush.

Compatible with all Dazzlepro Elite and Advanced model toothbrushes, the UV Sanitizing Charging Case is equipped with UV technology to effectively kill germs and bacteria in just 7 minutes while keeping your toothbrush fully charged for your next use.

The convenient USB charger makes it portable for use at home, office, or while traveling. Its slim, sleek design is ideal for counter tops, handbags, and carry-on luggage.

## Say Hello to a Beautiful Smile.



### Sanitizing Your Brush Heads

- Open the **UV Sanitizing Case Cover** and place your brush heads in their indicated positions. Make sure the bristles of brush head are facing the outside of the case.
- Close the **UV Sanitizing Case Cover**. The UV light will automatically activate to sanitize your brush heads and will shut-off automatically after 7 minutes at the end of the sanitization process. If the cover is opened accidentally during operation, the UV lamp will automatically turn off.

It is recommended to store your toothbrush and brush heads in the UV Sanitizing Charging Base to keep them charged and sanitized for each use.

### Safety & Precautions

- Basic safety precautions should always be followed when using electrical devices.
- If the device fails to function, do not attempt to dismantle, exchange or manipulate any part of it to avoid the risk of electric shock.
- This product is not intended for use by children under the age of three.


Contact Dazzlepro if you have any issues or safety concerns.



### USB Connecting and Charging

The UV Sanitizing Charging Case does not contain a battery. The USB cable must be connected to a power source in order to operate the unit.

- Insert the USB cable into the USB Power Port of the UV Sanitizing Charging Case.
- Plug the other end of the USB cable into a computer or USB power adaptor
- Place the toothbrush handle into the Charging Case so that the charging indicator lights up.
- For everyday use, keep the USB connected to a power source and store your toothbrush handle in the charging case to maintain its full power.

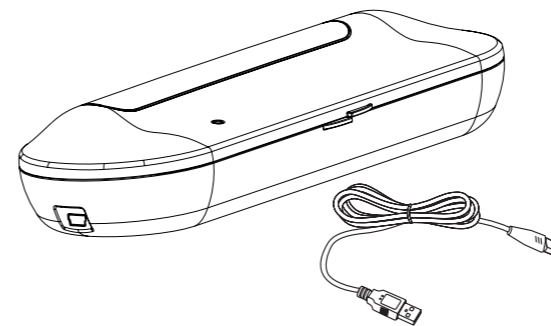

### Cleaning & Recommendations

- Dry your toothbrush handle thoroughly before placing it into the charging case.
- Place brush heads in their indicated positions. Make sure the bristles of brush head are facing the outside of the case.
- **Before cleaning your Charging Case, make sure to unplug the USB adaptor.** Use only a dry cloth to wipe the surface clean.
- **Store the Charging Case in a safe, dry place.**
- It is recommended to fully charge your toothbrush prior to first use. A full charge takes approximately 24 hours and will provide enough power for approximately 2 weeks of regular brushing.

**IMPORTANT:**  
Never place the UV Sanitizing Charging Case under running water or submerge it in water.