IMPORTANT SAFEGUARDS

READ ALL INSTRUCTIONS BEFORE USING

- DANGER

- DO NOT submerge the unit in water or any other liquids.
- WARNING

- Close supervision is necessary when this product is used by, on, or near children or individuals with certain disabilities.
 - Use this product only for its intended use as described in this manual. Do not use attachments unless recommended by the manufacturer.


SAVE THESE INSTRUCTIONS

MEDICAL WARNING

For your safety:

- Contact your dentist prior to use if you have had oral or periodontal surgery in the previous 12 months.
- When you start any new oral hygiene program, your gums may initially bleed slightly. Contact your dentist if excessive bleeding occurs or if bleeding continues after two weeks of use.
- Be sure that the water jet tip is held ¼ in. away from your teeth and gums. Damage to teeth and gums may occur if any other part of the unit touches them.
- This product is designed for oral cleaning only. DO NOT use the unit for any other purpose.
- As with any other new oral hygiene program, individuals with heart problems should check with their physician prior to using the Interplak® by Conair® Cordless Portable Water Flossing System, and should always exercise particular care with any oral hygiene device.
- The Interplak® by Conair® Cordless Portable Water Flossing System is a personal care unit and is intended for household use only.

SAVE THESE INSTRUCTIONS

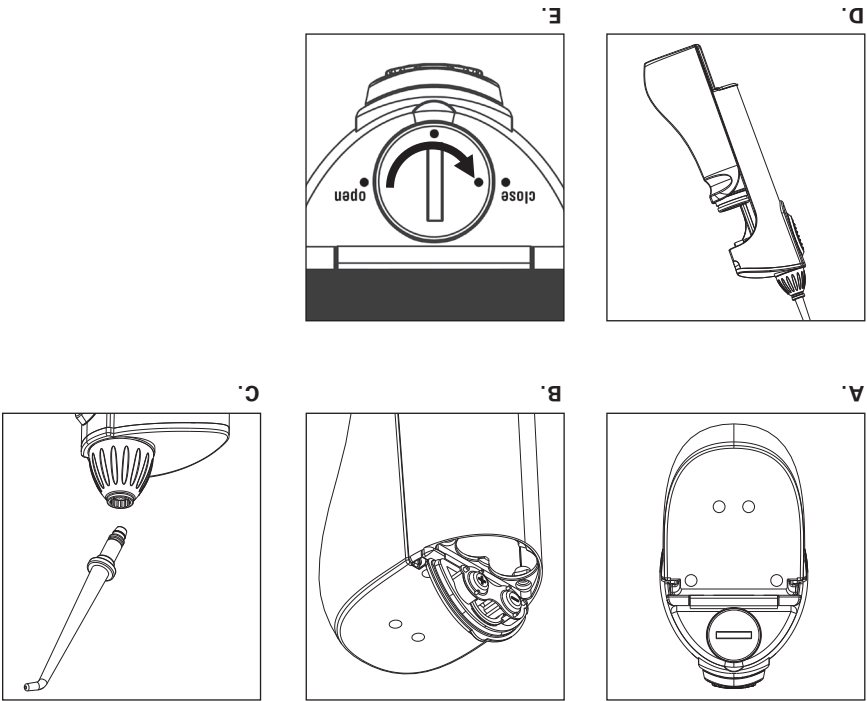NOTE: You can also use your favorite dental rinse or mouthwash. Thoroughly rinse reservoir with clean water after each use with mouthwash or dental rinse, and whenever unit has not been used for an extended period of time.

3. Begin with switch positioned to lowest setting (Lo) to allow gums to become accustomed to water jet cleaning. Before turning your unit on, place the jet tip in your mouth, keep your lips slightly closed, and **lean over the sink**. Let water flow from your mouth into the sink.

4. Adjust to higher setting (Hi), if desired, for more powerful waterflow. **NOTE:** If the pressure diminishes over time, replace batteries.

5. Turn handle to direct the water jet tip toward tooth surfaces and between teeth to rinse away plaque and trapped particles. The pulsing action of the water massages gums as it cleans the entire mouth.

NOTE: To prevent damage to delicate gum areas, do not point water jet tip directly at gums. Water jet tip should be held ¼ in. away from your teeth and gums.

1. Place water jet tip on unit by snapping tip into place (see C).

2. There is only one way to fill the water tank. Remove the water tank by sliding it down the unit (see D).

Initial Preparation:

Limited 1-Year Warranty (U.S. and Canada Only)

Conair will repair or replace (at our option) your unit free of charge for 12 months from the date of purchase if the appliance is defective in workmanship or materials.

To obtain service under this warranty, return the defective product to the service center listed below, together with your purchase receipt and \$3.00 for postage and handling.

California residents need only provide proof of purchase and should call 1-800-334-4031 for shipping instructions. In the absence of a purchase receipt, the warranty period shall be 12 months from the date of manufacture.

ANY IMPLIED WARRANTIES, OBLIGATIONS, OR LIABILITIES, INCLUDING BUT NOT LIMITED TO THE IMPLIED WARRANTY OF MERCHANTABILITY AND FITNESS FOR A PARTICULAR PURPOSE, SHALL BE LIMITED IN DURATION TO THE 12-MONTH DURATION OF THIS WRITTEN, LIMITED WARRANTY.

Some states do not allow limitations on how long an implied warranty lasts, so the above limitation may not apply to you.

IN NO EVENT SHALL CONAIR BE LIABLE FOR ANY SPECIAL, INCIDENTAL, OR CONSEQUENTIAL DAMAGES FOR BREACH OF THIS OR ANY OTHER WARRANTY, EXPRESS OR IMPLIED, WHATSOEVER. Some states do not allow the exclusion or limitation of special, incidental, or consequential damages, so the above limitation may not apply to you.

This warranty gives you specific legal rights, and you may also have other rights, which vary from state to state.

Please register this product at: www.conair.com/registration

Service Center:
Conair Corporation, Service Department, 7475 N. Glen Harbor Blvd., Glendale, AZ 85307

©2017 Conair Corporation, 7475 N. Glen Harbor Blvd., Glendale, AZ 85307
For more information regarding this and other Conair products, visit www.interplak.com or call 1-800-334-4031

71IN030525
IB-15283

Instruction Manual • Model WJL (All Versions)

For your safety and continued enjoyment of this product, always read the instruction book carefully before using.

Questions or Comments? 1-800-334-4031
www.interplak.com

Welcome to a New Level of At-Home Dental Care

The Interplak® by Conair® Cordless Portable Water Flossing System enhances oral hygiene and significantly contributes to healthy gums when used as part of a responsible oral health management program. It is particularly effective at cleaning around orthodontic appliances, fixed and removable bridge work, implants and splints, and it provides valuable gum massage with every use.

Remember: Oral irrigation is a supplement to, not a substitute for, your oral hygiene routine. Ask your dental health care provider to advise you on how the water jet can most effectively meet your specific oral hygiene needs.

Water jet tip

HI, Lo, Off switch

Battery door

Removable water tank

Before using your Interplak® by Conair® Cordless Portable Water Flossing System for the first time:

1. Turn battery door lock using a coin or flathead screwdriver and insert 2 AA alkaline batteries (see A). Ensure that the batteries are installed correctly with regard to polarity (+ and -) (see B). Replace door. Always purchase the correct size and grade of battery most suitable for intended use. Replace all batteries of a set at the same time.

To avoid water getting into the battery compartment, the door must be firmly closed and **securely locked**. Slip a coin or flathead screwdriver into the slot to lock and turn clockwise to tighten (see E).

NOTE: Use only alkaline batteries. Do not mix old and new batteries. Remove alkaline, standard (carbon-zinc), or rechargeable (nickel cadmium) batteries. Remove batteries from equipment that is not to be used for an extended period of time. When replacing batteries, make sure you dispose of the old batteries properly. Clean the battery contacts and also those of the device prior to battery installation. Remove used batteries promptly.

When replacing the batteries, make sure you dispose of the old batteries in accordance with applicable regulations. They should be disposed of in a location specifically designated for that purpose so they can be recycled safely and without posing any risk for the environment. Do not burn or bury them.

2. Prior to filling the water tank, make sure that the water jet tip is firmly placed onto the top of unit (see C).

3. Allow full reservoir of lukewarm water to pass through the tip before first use. (See Instructions for Use.)

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