

waterpik® waterflosser®


QUICK START GUIDE

WP-440/450/455/460/462/
463/464/465/466


Register your product and
receive over \$50 in savings.
waterpik.com/registration

Visit waterpik.com/howto450
to watch How To Video

1 Unit Prep

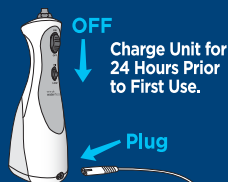

A. Turn power switch **OFF**, plug power cord into outlet, and insert the charging cord into the unit.

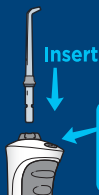


Caution:
Do not press tip
eject button while
unit is running.

B. Press tip firmly into unit.
To remove tip, push the tip
eject button.



C. Unplug charging cord.
Lift the flip top on reservoir
and fill with warm water.
Refill as needed.


2 Water Floss

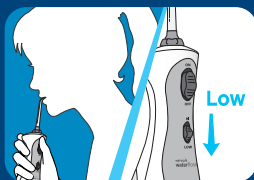

A. Start on low pressure and increase to your liking. With unit **OFF**, place tip in mouth. Turn power **ON**.



B. Lean over sink. Close lips enough to prevent splashing. Keep unit upright during use to ensure continuous water flow.

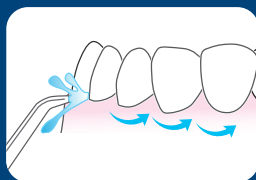

C. Aim water at the gumline at a 90 degree angle. Follow the gumline and pause briefly between teeth.

TO AVOID MESS: KEEP TIP IN MOUTH WHILE UNIT IS ON.

waterpik® waterflosser®

GUÍA RÁPIDA PARA COMENZAR WP-440/450/455/460/462/463/464/465/466

Registre su producto y reciba
más de USD 50 en descuentos.
waterpik.com/registration

Visite waterpik.com/howto450
para ver el video práctico

1 Preparación de la unidad



Off (Apagado)

Cargue la unidad
durante 24 horas
antes del primer
uso.

Enchufar

A. Colocar el interruptor de encendido en OFF (Apagado), enchufe el cable de suministro de energía en el tomacorriente e introduzca el cable de carga en la unidad.



Inserte

Precaución:

No oprima el botón
para eyectar la punta
mientras la unidad
está funcionando.


B. Presione con firmeza la punta dentro de la unidad. Para quitar la punta, presione el botón para eyectar la punta.



Llenar

C. Desenchufe el cable de carga. Levante la tapa abatible del depósito y llénelo con agua tibia. Vuelva a llenar si es necesario.

2 Cómo utilizar el Irrigador oral

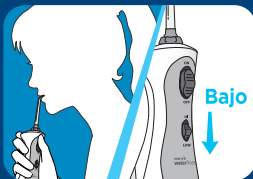


Bajo

A. Comience con presión baja e increméntela según su agrado. Con la unidad apagada, coloque la punta en la boca. **ENCIENDA** la unidad.



B. Inclínese sobre el lavabo. Cierre los labios lo suficiente para evitar salpicar. Mantenga la unidad en posición vertical durante el uso para asegurar el flujo continuo del agua.

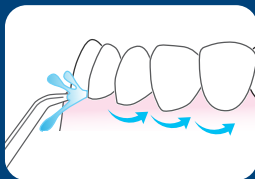

C. Dirija el agua hacia las encías en un ángulo de 90 grados. Siga las encías y haga pausas breves entre los dientes.

EVITE DERRAMES: MANTENGA LA PUNTA EN LA BOCA MIENTRAS LA UNIDAD ESTÁ ENCENDIDA.