
Let's
get
started

Meet your new Normatec Elite

Right leg

Left leg

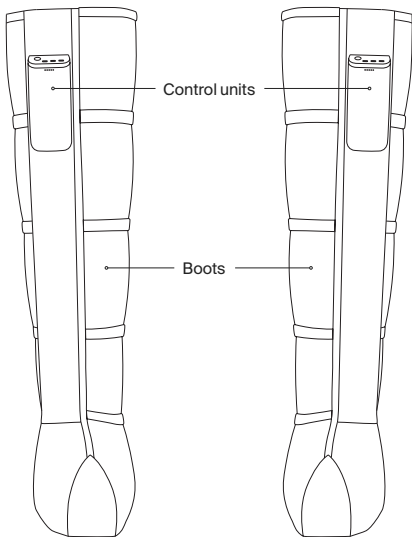Control unit

HyperSync™
status indicator

Bluetooth®
status indicator

Start/Stop button

Power (ON/OFF) button


Battery level indicators

ZoneBoost™ button

Time adjustment button

Pressure level
adjustment button

Power jack

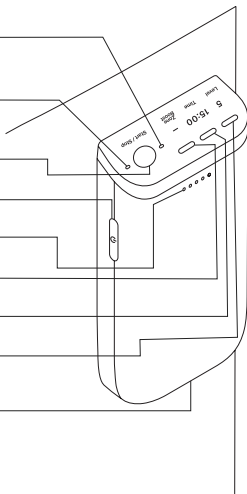Before you start


Download the Hyperice App

Get the most out of your Normatec Elite, or any Hyperice device with the Hyperice App. Remotely stop and start, adjust pressure level and time, and access advanced features and extra customization. Simply open the app to automatically connect to your device. Bluetooth® connection indicator will illuminate when connection is successful.

Register your device


Activate your warranty and ensure easy returns, repairs, or refunds at hyperice.com/register-product.

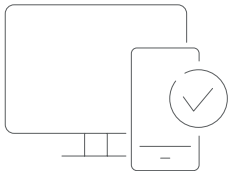Please read the Operating Instructions before first use.

Charge your new Normatec Elite


Plug in your Normatec Elite with the provided Hyperice power supply. Full charge is indicated when all five white battery status LEDs are illuminated and solid.

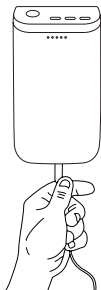Start your session


Putting on your Normatec Elite

The Normatec Elite can be worn over any comfortable clothing. Make sure your pockets are empty, and there are no bunches of fabric that will be uncomfortable during your session. Find a comfortable position sitting, reclining, or lying down, and put the Normatec Elite on be sure to zip up them completely. Never try to use the system with the zipper partially or totally unzipped.

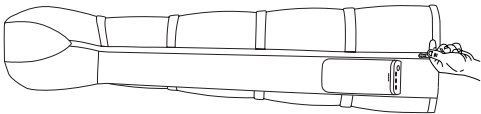
Powering up

Turn your Normatec Elite ON by pressing the power (ON/OFF) button. Hold for one second until display and battery level indicator illuminates.

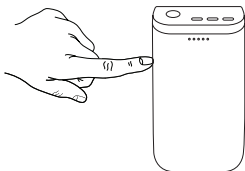
Adjusting pressure level

Choose your desired pressure by pressing once per level. The number on the display adjacent to the button indicates the current level.

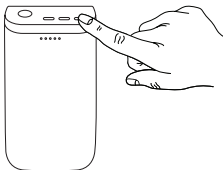
Adjusting treatment time

Choose your treatment time by pressing once per level (increments of 15 minutes). The numbers on the display adjacent to the button indicates the current setting.

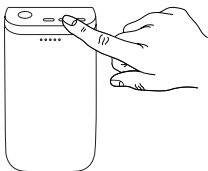
Adding ZoneBoost™

During your session, you can increase the intensity of a single zone with the ZoneBoost feature, when you want extra attention in a specific area. To boost a zone press the ZoneBoost button until the display shows the number of the zone you want to boost.

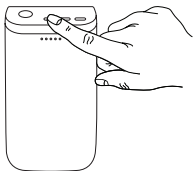Start your session

To start your session, press the Start/Stop button. To pause the session at any time, press the Start/Stop button. To restart your paused session, press the Start/Stop button again.

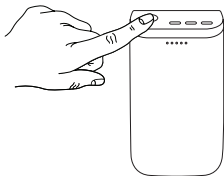
Ending your session

When your session has ended and treatment has stopped, allow the boots to deflate for 5-10 seconds before unzipping to prevent zipper damage.

Caring for your Normatec Elite

Make sure the power is OFF and the power supply is not attached.


Use a damp, clean cloth to gently wipe your device. Store in a clean, cool, dry place, out of direct sunlight, when not in use.

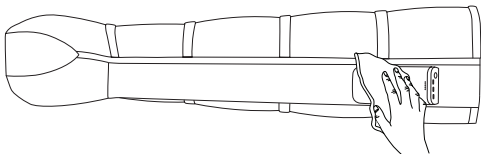Using HyperSync™


Unpair devices

Devices arrive paired, to unpair them press and hold the Start/Stop button on either device until the device screen reads “Unpaired!” and the HyperSync™ pairing indicator light is no longer illuminated.

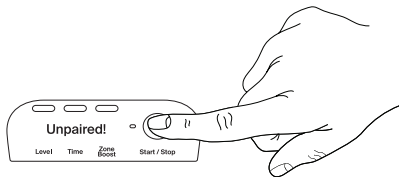Pair devices

Turn on both devices. Press and hold Start/Stop button on either control unit until the display screen says "Paired!" Press and hold the Start/Stop button on the other device until both devices say "Paired!" and the HyperSync™ pairing indicator light is illuminated.

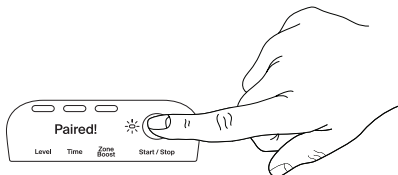We're here for you

HyperCare®

Get award-winning support from our HyperCare team – a group of specialists dedicated to your overall well-being and expert guidance on Hyperice products.

1.855.734.7224

hyperice.com

customersupport@hyperice.com

If outside of the United States, please visit hyperice.com/contact

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