CURAPROX

HYDROSONIC PRO

Manual

 SWISS PREMIUM ORAL CARE

CURAPROX HYDROSONIC PRO

HYDROSONIC PRO

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Congratulations

on buying the 'Hydrosonic pro' sonic toothbrush from Curaprox.

Curaprox is the oral health brand of Curaden AG, based in Kriens, Switzerland. Since 1972, Curaprox has been working together with experts from research, teaching and practice to develop gentle and effective tools for oral hygiene.

Hydrosonic pro

For professional teeth-cleaning at home. Especially suitable if you have braces, implants or gum problems and for keeping teeth and gums healthy every day.

Read this user manual carefully before you use the appliance and save it for future reference.

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Warning

If you've recently had surgery on your teeth or gums, consult your dentist before using the 'Hydrosonic pro'. If you experience excessive bleeding after using this device, please contact your dentist immediately. This toothbrush can be used by children and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved. Cleaning and user maintenance shall not be made by children without supervision. Children shall not play with the appliance. The 'Hydrosonic pro' complies with safety standards for electromagnetic devices. If you wear an implant or pacemaker and have questions about the device, please contact your doctor or the manufacturer of your implant or pacemaker.

Safety note

Warning: For the purposes of recharging the battery, only use the detachable supply unit provided with this appliance. This appliance contains batteries that are non-replaceable. Before using, please check that the voltage stated on the device complies with your local mains voltage. If the Hydrosonic charger or connecting cable is damaged, please do not use or attempt to repair it yourself. The power cord cannot be replaced. If the cord is damaged, the entire device must be replaced. All repairs must be handled by a dealer or our customer service. Keep the charger away from

water, and do not immerse in water or other liquids. Cleaning the charger: Disconnect from the mains before cleaning and reconnect only once the charger is completely dry. Never charge the device unattended.

Environmentally-friendly disposal

Do not dispose of this product with normal household waste at the end of its service life. Take it to a recycling centre or official collection point. This helps protect the environment. The integrated rechargeable battery contains substances that may endanger the environment. So please take the used handset to an official recycling centre or collection point.

Proper disposal of this product

Within the EU, this symbol indicates that this product must not be disposed of with household waste. Old appliances contain valuable recyclable materials. They should be recycled in order to prevent harm to the environment or human health through improper waste disposal. Therefore, please dispose of old devices through a proper collection system, or send the device back to where you purchased it for proper disposal. They will then ensure that the device is recycled accordingly.

Technical data for USB plug

Adapter input: 100-240 VAC, 50/60Hz, 0.2A max.
Adapter output: 5 VDC, 600mA

Technical data for charging station with USB cable

Input: 5 VDC, 600mA

Technical data for handset

- Battery: lithium-ion (3.7 VDC, 700mAh)
- Optimal battery charging time: 10 hours
- Operating time per battery charge: approx. 60 mins.
This is equivalent to 15 days if you brush your teeth for four minutes each day.
- Brushing movements per minute: 44,000, 64,000 to 84,000; the last two feature in three modes each
- Quadrant timer: 30 seconds
- Brushing timer: two minutes
- Automatic switch-off: after four minutes
- Accessories: 1× 'sensitive' brush head
1× 'power' brush head
1× 'single' brush head
1× charger with USB cable
1× USB plug
1× travel case
1× instruction manual

Troubleshooting

Problem	Possible cause	Solution
• The device isn't working	• The device is not turned on	• Place the handset securely on the charger • Charge the 'Hydrosonic pro' completely before first use. The charging takes about 10 hours.
• Brushing performance is getting weaker	• The brush head is dirty or worn	• Clean or change the brush head
• You feel a slight tingling when using the sonic toothbrush	• If you are using a sonic toothbrush for the first time, this is a normal sensation	• You will become accustomed to the sensation after using the sonic toothbrush for a while


HYDROSONIC PRO

The basics	34
Charging status	35
The seven modes	36
Using the brush and getting it going	37
Brushing with braces: Instructions	38
Brushing after braces: Instructions	39
Cleaning the gum line: Instructions	43
The single attachment: More instructions	44
Background: Mechanical and hydrodynamic	46
Background: Smart brushes for clever care	47
Reorder: Brush heads	50
Customer Service	51

30 DAY
**MONEY
BACK**
GUARANTEE

Contents

- 1× 'Hydrosonic pro'
- 1× 'sensitive' brush head
- 1× 'power' brush head
- 1× 'single' brush head
- 1× charger with USB cable
- 1× USB plug
- 1× travel case
- 1× instruction manual

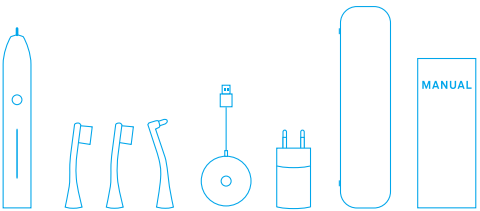
THE BASICS

Choose your brush head

- 'sensitive' and 'power' for cleaning healthy teeth and gums, as well as for gum problems, implants and orthodontic devices.
- 'single' for contact points with orthodontic devices, from brackets and lingual brackets to retainers; and for caring for the gum line.

The brush heads are easy to attach and remove.

Caution: Do not bend.

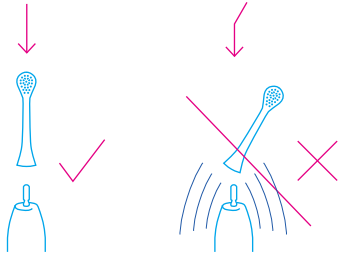Switch on and off Press button briefly.

Press briefly:
on/off

Mode selection


Press and hold button; release at desired mode. When you turn it on, the 'Hydrosonic pro' will start on the last selected mode. Bristle movement will automatically and immediately slow down when too much pressure is used.

Press and hold:
Mode selection

CHARGING STATUS: ALWAYS AT FULL POWER


- Charge the 'Hydrosonic pro' completely before first use.
- Optimal battery charging time: 10 hours.
- To charge, place the 'Hydrosonic pro' into the charger, which should be connected to the mains: plug the USB plug into the charger's USB connection, and plug into the mains. You can also charge the device using the USB port on your laptop or computer.
- You can see the charge status when you remove the 'Hydrosonic pro' from the charger or switch it off.
- When all seven LEDs on the display light up, the battery is fully charged: the more the LEDs light up, the more charged the battery.
- If the bottom LED is red instead of blue, the battery is less than 20 % charged.
- Once fully charged, the battery will last a good 60 minutes, which equates to two weeks of brushing for four minutes each day.
- For extra-long life: Charge your toothbrush only when the battery status falls below 20 % (the red LED will indicate this).
- Always at full power: cleaning performance does not change as the battery runs out. The 'Hydrosonic pro' cleans just as well on low battery as when it's fully charged.

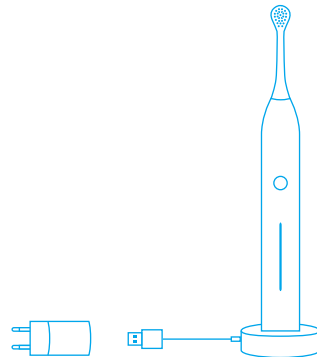Charging: Use the charger or the USB port on your laptop/computer.

THE SEVEN MODES: TRY THEM OUT

Switch on then hold the button: the 'Hydrosonic pro' will run through each of its seven modes.

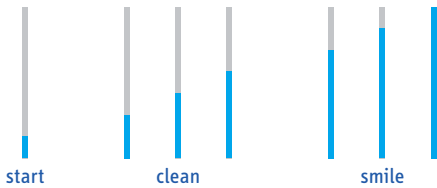One «start» mode

We recommend using the «start» mode for the special cleaning of contact points on orthodontic devices and prostheses, around implants and for caring for the gum line. Use the 'single' attachment. This mode features over 44,000 brushing movements per minute. The bottom LED is lit up.

Three «clean» modes

Continue to «clean» mode, which has three different intensity modes each with 64,000 brushing movements per minute; a new LED will light up for each mode. Recommendation: use the 'sensitive' or 'power' attachment.

Three «smile» modes

Finally, we have three different «smile» modes each with 84,000 movements per minute. An LED will light up for each mode. After that, the modes start over again. Recommendation: Use the 'sensitive' or 'power' attachment.

Adjustable speed

Select your favourite mode by releasing the button when you have arrived at the mode that feels best for you. It will be saved for the next time you brush.

Interruptions

Three short, consecutive pauses of 0.2 seconds each will indicate that 30 seconds have elapsed. A long interruption means that two minutes have passed. After four minutes, the 'Hydrosonic pro' turns off automatically.

Tip: Every mode delivers a great clean. With the exception of the 'start' mode for cleaning orthodontic devices and implants and caring for the gum line, all modes are strong enough to provide not just a mechanical clean, but also a hydrodynamic one.

USING THE BRUSH: THE BEST METHOD

- The brush heads are high-tech, and come with the finest Curen® filaments. Due to the 84,000 brushing movements per minute, there is an extraordinary amount of stress on the 'Hydrosonic pro'.
- So, you should change your brush head at least every three months.
- After using, rinse the brush head and bristles under running water.
- If you wear orthodontic equipment such as brackets or braces, the amount of wear is even higher. In this case, it is best to change your brush head after seven or eight weeks. It is worth it.

The 'Hydrosonic pro' looks great. Here's how to keep it looking that way. After use, simply hold it under running water and dry with a soft cloth. – If the handle is very dirty, wash it with a bit of soap and then dry it. Storage: avoid contact with water when storing – it's best to store it in the bathroom cupboard or on a dry surface. This will keep your 'Hydrosonic pro' looking great for a long time.

GETTING IT GOING: WITH FEELING

- Choose your favourite level from 'clean' to 'smile' while holding the brush to your tooth – do it by feel. This ensures brushing will be especially pleasant.
- And yes, try both brush heads, 'sensitive' and 'power'. Which feels best?
- From a cleaning perspective, it doesn't matter which of the three 'clean' or 'smile' levels you use.


Note: The 'clean' and 'smile' modes are suitable for cleaning healthy teeth and gums, as well as for gum problems, implants and orthodontic devices – we recommend using the 'single' brush head only in the 'start' mode.

BRUSHING WITH BRACES: INSTRUCTIONS

Brackets affixed on the teeth create additional niches to clean. Taking good care of these niches is a deciding factor in oral health. Whether you have brackets, lingual brackets or retainers – thanks to the 'Hydrosonic pro', brushing is easy and your teeth stay healthy.


Basic cleaning

Using the 'sensitive' or 'power' brush head, go from tooth to tooth between your bracket and gum line, so above or below the bracket only, in one of the «clean» to «smile» levels. This means you're also cleaning your gum line at the same time. Finally, clean the other side of your bracket, moving from tooth to tooth again.

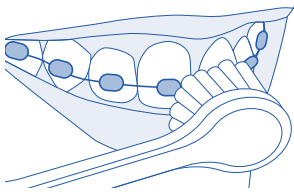
Brackets

How to clean contact points: Using the 'single' attachment, apply minimal movements to go around the brackets, nice and slowly, with almost no pressure, in the 'start' mode.

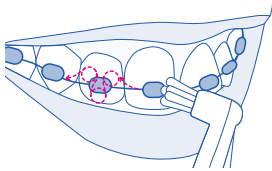Lingual brackets and retainers


The same goes for brackets secured to the backs of your teeth. Using the 'single' attachment, go nice and slowly along the lingual brackets or retainer with almost no pressure, in the 'start' mode. Thanks to Curacurve®, these areas on the backs of your teeth are easy to reach.

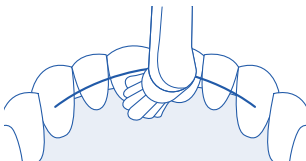BRUSHING AFTER BRACES: INSTRUCTIONS

And after braces? Carry on cleaning thoroughly. You have everything you need: the 'Hydrosonic pro', the 'sensitive' and 'power' brush heads and the 'single' attachment.




1. Moisten brush head and apply a little toothpaste.
2. Put the brush in your mouth and place it at a 45-degree angle so that the brush is lying gently half on the teeth and half on the gums.

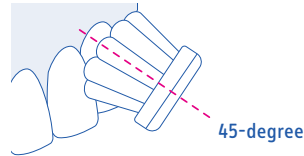
3. Clean from tooth to tooth: Stay on each tooth at a 45-degree angle for two to three seconds. Do not apply any pressure.

Tilting motion



Tilt the 'Hydrosonic pro' once or twice, as if trying to insert the bristles between the teeth. Do this very slowly and without applying pressure.

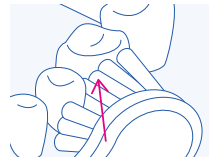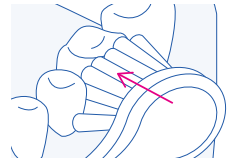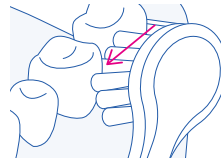Remain on each tooth for two or three seconds while tilting the brush once or twice to thoroughly clean away build-up.


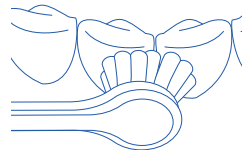
4. After 30 seconds, three short, consecutive pauses of 0.2 seconds each will indicate that it's time to move to a different quadrant of your mouth.
5. Repeat the cleaning process for all quadrants. Make sure that you always use the same order when brushing so that you don't skip a tooth.

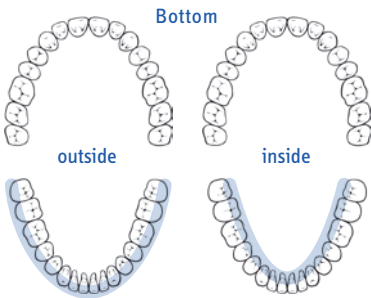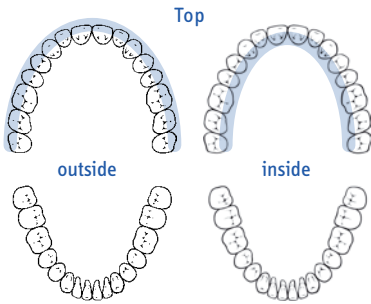6. A long interruption indicates that two minutes have passed. After four minutes, the 'Hydrosonic pro' turns off automatically.
7. You clean permanent implants with the 'Hydrosonic pro' just as you would clean natural teeth. You don't need to worry if you have an implant. We highly recommend caring for implants with interdental and/or implant brushes, because the gum line around the implant should be cleaned with particular care. Please consult your dentist.

For hybrid implants: the Hydrosonic is ideal for cleaning bars and buttons.

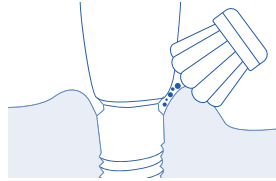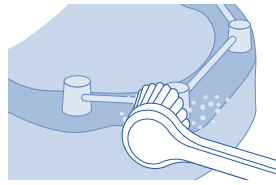8. It's especially important that children and young people also clean chewing surfaces that may have deep trenches. Just as before and using no pressure, move the 'Hydrosonic pro' from tooth to tooth, from chewing surface to chewing surface, and let the 'Hydrosonic pro' do its work for two to three seconds each time.



9. For fillings, implants, crowns and veneers: Simply use the process described in these instructions.
10. According to dental recommendations, you should brush for at least two minutes. Extending this cleaning time will improve oral health. The toothbrush may be switched off at any time.

Important


Do not exert any pressure on the toothbrush. To ensure optimal cleaning, place the 'Hydrosonic pro' gently onto the surface of the tooth and clean tooth by tooth. This approach will ensure optimal prophylaxis. If your gums become sensitive after the second or third use, you are not using the 'Hydrosonic pro' properly. Please read these instructions again carefully or consult your dental hygienist, dentist or orthodontist for advice.

Special cases

- a. Upper and lower front teeth: To clean the inside part of the lower and upper front teeth, hold the 'Hydrosonic pro' vertically and tilt it to the right and to the left a few times before moving on to the next tooth.

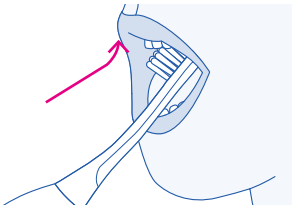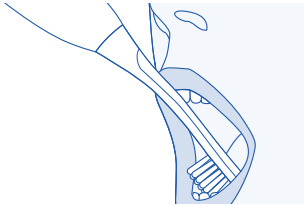- b. Upper and lower molars: Place the 'Hydrosonic pro' at an angle behind the molar and linger for two to three seconds.

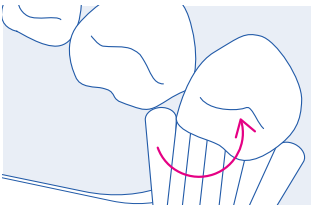CLEANING THE GUM LINE: INSTRUCTIONS


The 'single' attachment should really be called the «multipurpose» attachment, as it has so many different uses.

Cleaning the gum line


A lot of bacteria hides in the narrow furrow between the teeth and gums. The 'single' brush head allows you to clean the gum line intensively yet gently:

1. Place the 'single' brush on one side of the tooth, perpendicular to the tooth's surface.


This doesn't require toothpaste.

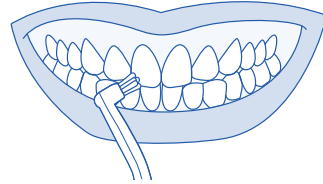2. The outermost bristles should be able to gently and sensitively enter the gum line a little almost without any pressure at all.

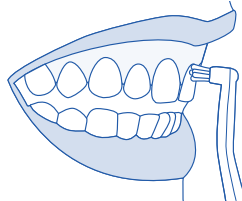
3. Now, move the 'single' brush along the gum line without applying any pressure, moving from one side of the tooth to the other. Always hold the 'single' brush head perpendicular to the surface of the tooth.

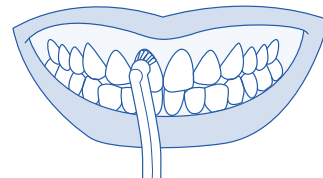HOW TO USE

4. Proceed in the same way for the back side of the tooth. Place the 'single' brush head at one end of the tooth. Touch the gum line gently and clean the tooth using the same principle you used on the front.

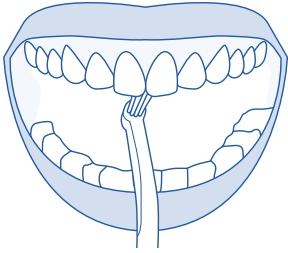
Now move on to the next tooth. Then use either the toothbrush attachment or the 'single' brush head to clean the chewing surfaces.

THE SINGLE ATTACHMENT: MORE INSTRUCTIONS


Awkwardly positioned teeth

For awkwardly positioned teeth, move along the niches and overlaps with other teeth slowly, almost without any pressure at all.

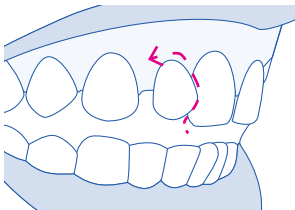Extensive tartar build-up

People with especially extensive tartar build-up should move the 'single' brush along the gum line, applying almost no pressure at all, in the 'start' level.

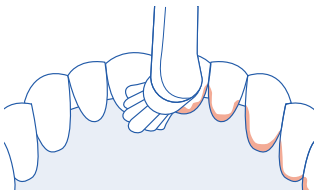
Implants

To ensure a long life, implants need good care.


Clean critical areas particularly well and almost without any pressure: the gum line around the implant.

The single brush head is the best tool for this.

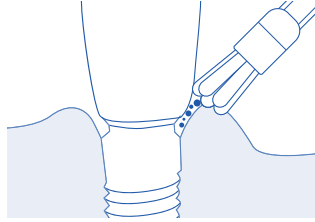Polishing

Clean each tooth surface with the 'single' brush head, almost without any pressure, and use small circular movements in the 'start' level. Use your tongue to feel how smooth and clean your teeth are.

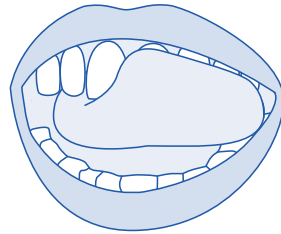BACKGROUND: MECHANICAL AND HYDRODYNAMIC

Effective and gentle

With over 60,000 brushing movements per minute, the 'clean' and 'smile' modes let you clean more effectively and gently than with a manual toothbrush: It's almost impossible to damage your teeth or gums.

Your teeth and in between

Thanks to sonic technology, the 'Hydrosonic pro' provides an unbelievably thorough clean, not just on teeth and the gum line, but also to an extent in the spaces between your teeth.

Gum line, pockets

The 'Hydrosonic pro' removes millions of biofilm bacteria from the gum line and possible periodontal pockets – thanks to its sonic technology and hydrodynamics.

Mechanical cleaning

The speed of the bristle movements ensures excellent mechanical cleaning performance.

Hydrodynamics

This high speed leads to turbulence in the mixture of saliva and toothpaste in your mouth. This turbulence removes lots of biofilm bacteria where the bristles cannot reach, like the spaces between your teeth.


Together with well-designed brush heads, the 'Hydrosonic pro' fulfils the highest demands in terms of mechanical and hydrodynamic cleaning.

BACKGROUND: SMART BRUSHES FOR CLEVER CARE

How good a sonic toothbrush is depends on how powerful it is – and even more important are the brush heads: both the arrangement and quality of the bristles.

Filaments


We have been using Curen® filaments rather than standard nylon bristles from the very beginning. The benefit: Curen® hardly absorbs any water, so performance remains constant. And especially important: Curen® is incredibly soft and feels wonderful – even at 84,000 brushing movements per minute.

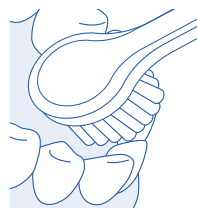More on page 39

Arrangement


Professor DMD Ulrich P. Saxer recognised the benefits of the drop shape: it ensures each tooth is cleaned individually, resulting in a high-quality clean. The bristles are also arranged at different heights so that teeth are almost enveloped.

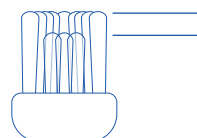Long on the outside, short on the inside: this filament arrangement makes sure different surfaces of each tooth are reached for optimal mechanical and hydrodynamic cleaning.


Movement

The circular motion of oscillating brush heads can damage the gum line if too much pressure is applied, the bristles are too hard, or if other improper techniques are used. That's why our oral health professionals opted for a simple vertical to and fro. This has three benefits: It's gentle on the gum line. It provides especially good care. And it provides an additional cleaning approach alongside the mechanical – the hydrodynamics. The filaments move so quickly that the turbulence created moves the mixture of saliva, water and toothpaste in the mouth. This means that areas the filaments can't reach, such as in between the teeth, are also cleaned.

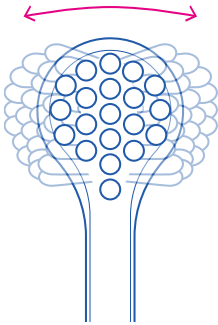The special case

Constructed by oral health expert Jiří Sedelmayer as the manual CS 'single' toothbrush, our 'single' attachment is now also available for the sonic toothbrush. The 'single' brush is extremely beneficial as a sonic toothbrush attachment: it lets you easily clean niches and furrows, such as along the gum line and around brackets, to perfection.

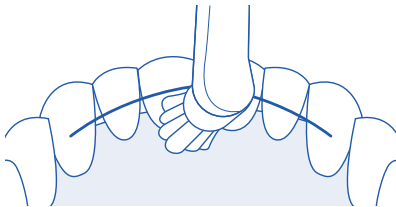Ergonomics

Yes, our bristle attachments are slightly bent. This curve of 10 or 15 degrees makes a difference: Ergonomics are considerably improved, difficult areas such as behind molars or the backs of the front teeth are easily reached – thanks to Curacurve®, the perfect curve.


10 degrees

15 degrees

0

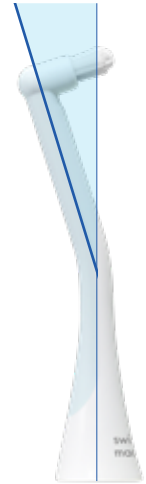0

Comfort

How can we make your brushing experience even more pleasant? Sonic toothbrushes in particular can feel unpleasant when the back of the brush head touches a tooth. Our solution: We have rubberised the back of the brush in a way that does not hinder the brush head's performance. Pleasant, simple, clean.

Our innovations

We are proud of our innovative Swiss brush head technology – and that we produce our brush heads in Switzerland. And we love it when smart users appreciate our efforts. We are convinced that our brush heads will help you achieve optimal oral health easily and comfortably.

MADE IN SWITZERLAND

REORDER: BRUSH HEADS

You can buy brush heads in your pharmacy, chemist's or dental practice or on curaprox.com.

power

Powerful yet gentle. With our soft Curen® filaments and Curacurve®, the perfect curve. For normal gums, or just because you like it. Extra small drop-shaped head for highly precise cleaning, tooth by tooth. Supports the hydrodynamic effect.

sensitive

Especially gentle yet highly efficient. With extra soft Curen® filaments and Curacurve®, the perfect curve. For sensitive gums, periodontitis, or just because you like it. Extra small drop-shaped head for highly precise cleaning, tooth by tooth. Supports the hydrodynamic effect.

single

Extra fine. Cleans in areas where bacteria are particularly well hidden: in the niches between teeth, brackets and wires. And the gum line, too. Almost essential for caring for implants and for reaching poorly positioned teeth. And for polishing tooth surfaces.

CUSTOMER SERVICE

Guarantee terms

Curaden AG provides guarantee for material, function and manufacturing defects for two years from the date of purchase. This may be carried out by way of repair work or a replacement device. In order for a guarantee claim to be made, the original purchase receipt must be submitted. Keep your receipt as proof of purchase. The guarantee shall expire if the device is subject to interference or changes to its original state.

Guarantee exclusions

In the following cases, Curaden AG shall not provide guarantee services: Wear and tear of the device or parts of the device due to normal use, or if the device was not used, not properly maintained, or if the safety instructions provided with the device were not observed. Repairs may only be carried out by service partners authorised by Curaden AG.

Service

In order to make a guarantee claim, send the device and original receipt to your country's official Curaprox representative. If you do not know who your country's representative is, visit our website www.curaprox.com or contact us directly in Switzerland: ask@curaprox.com.

- Please discuss any questions about dental or orthodontic care with your dentist. You can also contact us, especially if you have technical problems or if something is not working: ask@curaprox.com.
- Our FAQs are often helpful: www.curaprox.com/faq. Here, you will also find out how to exchange your 'Hydrosonic pro' if necessary.
- The 'Hydrosonic pro' isn't for you? You can return it for a full refund within 30 days of the purchase date.
- The guarantee expires after two years as standard. Register your device to get another one year guarantee for free: www.curaprox.com/warranty

These brush heads were developed together with Professor DMD Ulrich P. Saxer and oral health specialist Jiří Sedelmayer.

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