EN

With your inspiring rating, COSTWAY will be more consistent to offer you EASY SHOPPING EXPERIENCE, GOOD PRODUCTS and EFFICIENT SERVICE!

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
US office: Fontana **UK** office: Ipswich

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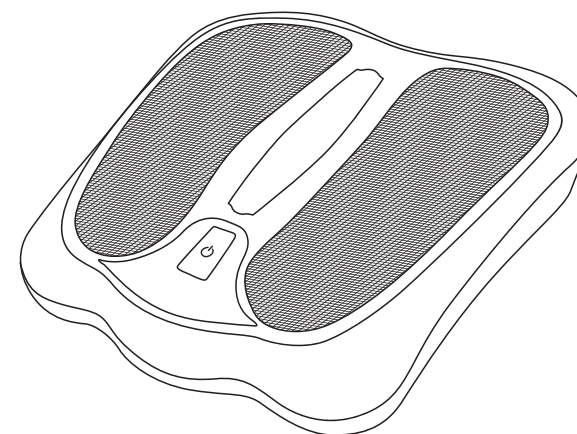To Users

Dear user:

We would like to convey in this letter our heartfelt thanks to you for your buying our production named far-infrared &kneading foot massager. Please read its instruction carefully before using the massager, especially safety attentions, so that you can use it correctly.

Besides, please keep the instruction after you read it in case you need refer to it someday.

Thanks for your cooperation!

Safety Attentions

1. Persons who are not allowed to use the massager

 WARNINGS	• Patients whose foot has skin disease. • Pregnant women, patients with heart disease.
--	---

2. Something you have to know when you are using the massager.

 ATTENTIONS	• Children are not allowed to use the massager unless with the adults' help. • Persons who are not sensitive to heat can't use the massager for toolong every time. • Don't stuff anything to the open parts. • Don't pull the massager on the ground in case it may cause damage to massager. • Don't put your center of gravity on one feet. • The massager is intended for persons'life and it can't be used in other ways.
--	---

 FORBIDDEN	• Stop using the massager if you find it damaged. • Don't stand on the massager or put something heavy on it. • Stop using the massager if it is abnormal.
 NO DEMOUNTING	• Don't take apart the massager by yourself, or it may break down. • Please contact with our after-sales service department if you have any questions.
 PULL OUT	• Don't pull out the attaching plug with wet hands. • Pull out the attaching plug after you using it. • Pull out the attaching plug when there is a power cut.

3. Concerning the attaching plug and power line


 ATTENTIONS	• Don't use the power supply whose voltage is higher than rated voltage, or it may cause fire. • Pull out the attaching plug instead of power lines. • Do use the dedicated power lines to replace the damaged one.
--	---

4. Operating environment


 ATTENTIONS	• Don't put the massager in a wet or dusty place. • Don't put the massager near somewhere in a high temperature. • Don't use the massager when it thunders.
--	---

Functions

1. Sketch map

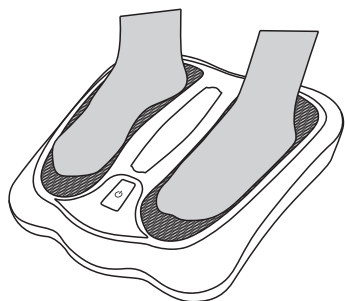2. The name of parts

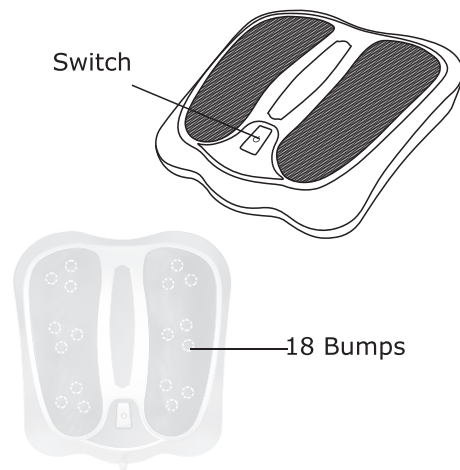3. Features

- Deep kneading and shiatsu massage
It can stimulate foot reflection zone to relieve the pain and relax the channels and collaterals.
- New massage head with 18 rotatable bumps.
You can enjoy the massage as if a professional massagist does it.
- Thermal therapy with far infrared
It helps to promote cell metabolism.

4. Press the switch twice and you can enjoy the kneading massage with heating. (The red light will shine.)
5. Pressing the switch three times means turning off the massager.

! Attentions

- Don't use the massager over 20 minutes. (It has no timing devices.)
- You'd better use the massager with socks on foot.
- Heating cannot be used alone.

Maintenance

1. Cleaning

- Before cleaning, pull out the attaching plug and wait until it cools down.
- Please clean the shell of the massager with a dry rag.
- Don't immerse the massager in any liquid.
- Don't use brushes, petrol, cleanser or diluent during cleaning.

2. Storage

- Put the massager at a safe and airy place if you won't use it for a long time.
- Don't hang the massager with wires.

Usages

Operation

1. Insert the attaching plug before using. (*alternating voltage : 220-240 V
Frequency: 50/60Hz*)
2. Sit down and put your foot on the massager.
3. Press the switch once and you can enjoy the kneading massage. (*The blue light will shine.*)

Trouble Shooting

If you meet problems when you are using the massager, please try to resolve it according to the tips below.

Abnormal situations	Possible reasons	Solutions
It doesn't work.	1. It can't energize. 2. You don't turn on the switch.	1. Insert the attaching plug again. 2. Turn on the switch.
Stop working suddenly.	The electricity is cut off.	Check whether the power intertace is in order.

If you still cannot resolve the problem, please pull out the attaching plug and contact with our after-sales service department.

Accessories

- ▶ 1x Foot Massager
- ▶ 1xUser's Manual

Technical Parameters

Model: HW66863DE/Z303

Power Rating: 31W

Voltage Rating: 220-240 V (AC)

Frequency Rating: 50/60Hz

Keep it in case you need it.

It may cause some pollution if you don't deal with the waste correctly.

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