

NuvōMed®

SLM-6/0110

**Electronic Shiatsu Heating Massager
for Neck, Shoulders and Back**



Instruction Manual

Please read all instructions carefully and retain for future use

Getting Started

- Remove the product from the box.
- Please read all instructions and retain all packaging until you have checked and counted all of the parts.


Package Contents

- Shiatsu Heating Massager
- 12V Power Adapter
- 12V Car Plug
- Instruction Manual

Key Features

- Relieve muscle tension and soreness to remove fatigue
- Make your body feel more relaxed
- Total 8 Massage Nodes (4 Nodes on Each Side)
- Dual use: Automobiles and Home
- Clockwise and Counter-clockwise massage
- One Key Operation
- 360° Deep Kneading

Description of Parts

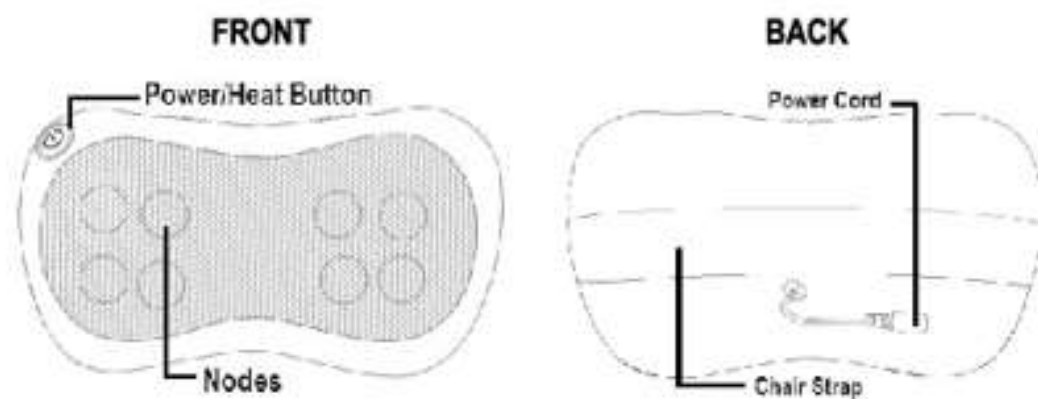

How to Use

- **Using 12V Power Adapter:** Plug the AC Adapter into the Power Cord end. Then, plug the input side into an electric outlet.
- **Using 12V Car Plug:** Plug the Car Plug into the Power Cord end. Then, plug the cigarette plug into 12V cigarette lighter socket inside your vehicle.
- The massager will automatically shut off after 15 minutes of continuous operation.
- Allow the massager to cool down for 30 minutes before using again. This will prolong the life of the massager.

Mode: On with Heat → Reverse Rotary Direction → Massage without Heat → Off

- Press the button once to turn the massager and heat on.
- Press again to reverse the rotary direction of massager.
- Third press to turn the heat off. Massage nodes keep working in a forward/reverse direction
- Forth press to turn the massager off.
- Unplug the massager when not in use.

How to Store

- When not using the massager, store it back in the box and place it in a dry, well-ventilated, cool area. Keep the product clean and away from moisture.
- **DO NOT** twist wires.
- Avoid sharp objects to prevent scratching surface of product.
- Keep out of direct sunlight and high temperatures.

How to Clean

- Unplug power cord before cleaning.
- Use a dry cloth only to spot clean the surface of the product.
- Avoid direct contact with any liquid or cleaning solutions.
- **DO NOT** use corrosive cleaners such as benzene or bleach products.
- **DO NOT** immerse the product in water.

Product Specifications

Product Size: 12" (L) x 7.5" (W) x 4" (H)
Voltage: 12V 2000mA

Product Weight: About 35.1 Oz
Rated Power: 24W
Rated time: 15 minutes

General Information

- Keep cord away from heated surfaces.
- Never block the air openings of the product or place it on a surface, such as bed or couch, where the air openings may be blocked. Keep the air openings free of lint/hair.
- Never use while sleeping.
- Never drop or insert objects into any opening.
- **DO NOT** operate where aerosol products are being used/administered.
- **DO NOT** use on/near eyes or other highly sensitive areas.
- **DO NOT** operate under a blanket/pillow. Excessive heating can occur and cause fire, electric shock/injury.
- **DO NOT** carry this product by its supply/power cord, or use the cord as handle.
- **DO NOT** use it directly on one concentrated area of the body. It is meant to be used in a continuous motion, constantly changing the massage area.
- To disconnect, turn off the unit then remove the plug.
- **DO NOT** perform any servicing on the product.
- Use this product only for its intended use as described in the manual. Do not use attachments not recommended by the manufacturer.

Safety Precautions

- Consult your doctor before using this product.
- Individuals with pacemakers should not use this item.
- This unit is not intended for use by persons (including children) with reduced physical, sensory or mental capabilities or lack of experience or knowledge. They should be given supervision and instruction in the use of the unit by a person responsible for their safety.
- Never use the product on open skin eruptions, swollen, or inflamed areas.
- DO NOT operate after taking painkillers, sedatives or alcohol. Consult your doctor if you are under these medications.
- DO not continuously operate for 15 minutes or longer. It is recommended to turn the product off and allow it to rest (cool) for 30 minutes before re-using. This will prolong the life of the product.
- The product is designed as a non-professional massager to soothe tired and aching muscles. It is not intended to be a substitute for medical attention. If your ailments and symptoms persist or worsen after use, contact your doctor immediately.

WARNING: Any individual who is pregnant, has a pacemaker, suffers from diabetes, phlebitis, and /or thrombosis, or is at an increased risk of developing blood clots due to recent surgery, should consult their physician before using a massaging device designed for home use.

FCC Instructions

Caution: Changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the product.

NOTE: This product has been tested and found to comply with the limits for a Class B digital device, pursuant to Part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This product generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If the product does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the product and receiver.
- Connect the product into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

The device complies with Part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) this device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.